



TOKYO JAPAN

English User's Manual come with this product.

Please download the other language's User's Manual version from this link (<https://www.ya-man.co.jp/en/support/documents/>) or scan with QR code with your camera phone.

Please carefully read User's Manual to ensure safe and proper use of the device.



TO ACCESS SCAN WITH
YOUR CAMERA PHONE

User's Manual

Home-care Beauty Device

Blue Green Mask Lift

YJMF4

Introduction

Thank you for purchasing a YA-MAN product.

- Do not use the device before reading the Safety Precautions (pages 5-11).
- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.
- This appliance contains batteries that are non-replaceable.

Steps for Use


1

Care for your skin
after ensuring
that your hair
does not fall onto
your face



2

Wear the mask
on your face.

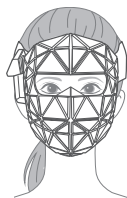

***Examples of incorrect wearing**



There is a space
between the mask
and your forehead


There is a space
between the mask
and your cheek or chin

There are gaps, and the electrodes are not
touching

To wear the mask properly

1

How to apply the gel

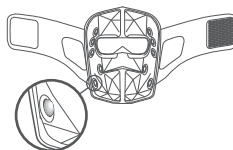

*Use about one pearl-sized drop of
the included skincare product for
each electrode.

2


How to wear the mask



- Attach by pulling up from under the chin
to the ○ section.

- Attach in a way there are no large gaps
between the mask and your face.


If the cheek or chin sections (● parts) lift up, the
electrodes near the jaw muscles will not be in contact
with the skin

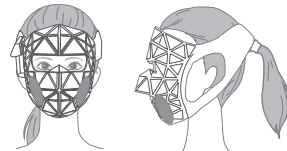

- It is fine even if the mask fits closely
to your mouth.

If you have a small face

Please use a headband or adjust with a hair tie.

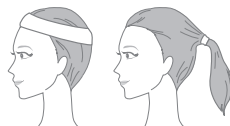


Steps for Use

3

Turn on the device

Press and hold the Power/LED Mode button for about two seconds.


Power/LED Mode button

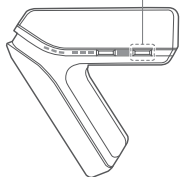

4

Select the LED Mode and the EMS level

Power/LED Mode button
EMS level button

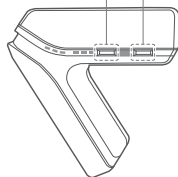

5

Turn off the device

The device will automatically turn off about 10 minutes after it has been turned on.

Precautions

- **Handle with care after use.**

After use, wipe off any gel remaining on the electrodes and allow them to dry completely.

Do not leave the mask in places where it is not easily visible. If left in such places, you may not notice it and accidentally step on it, which may cause device malfunction.

- **Do not apply excessive force to the mask.**

Do not apply excessive force to the mask, such as stepping on it, hitting it, or forcibly bending it.

Failure to comply may cause device malfunction.

- **Do not place the mask in ways that put a strain on it.**

The product will last longer if you store it folded with the controller on top.

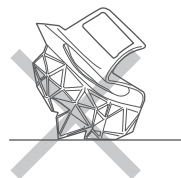
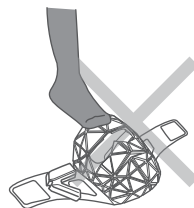


Table of Contents

Safety Precautions ————— p.5

Names and Functions of Components — p.12

Before Use ————— p.14

How to Use ————— p.16

After Treatment ————— p.20

Additional Information ————— p.21

● Troubleshooting ————— p.21

● Frequently Asked Questions ————— p.23

● Product Specifications ————— p.25

After-sales Service ————— p.26

Safety Precautions

- Read these Safety Precautions carefully before using the device in order to use it correctly.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- Do not use the device in a defective state, as doing so could cause a physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or CAUTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

Classifications of Precautions



DANGER

Improper use of this device poses a severe and likely risk of death or serious injury.



WARNING

Improper use of this device may pose a risk of death or serious injury.



CAUTION

Improper use of this device may result in injury or property damage*.

Symbols



This symbol shows you what you must not do.



This symbol shows you what you must do.



This symbol shows you what you must be careful with.

*'Property damage' refers to damage affecting dwellings, household belongings, or domestic animals and pets.

Danger (Use)

- **The following individuals should not use the device:**

Failure to comply may cause a malfunction in the medical electronic device.

◆ Individuals using the following types of electronic medical devices:

- pacemakers or other implanted electronic medical devices;
- artificial heart-lung machines, or other electronic life-support medical devices; and
- electrocardiographs or other electronic medical devices worn on the body.



Prohibition

WARNING (Use)

- **Do not use under the following conditions:**

Failure to comply may cause accidents or other problems.

- If you have existing or suspected heart disease or related disorders
- Infants or those who are unable to express their intentions on their own
- If you are unable to sense temperature or stimulation
- If you have an autoimmune disease

Consult a physician prior to use under the following conditions:

- If you have an acute disease
- If you have an infectious disease
- If you have a malignant tumor
- If you have a febrile disease
- If you have high blood pressure or a blood disease
- If you have a fever
- If you have sensory impairment due to severe peripheral circulatory dysfunction (e.g., diabetes)
- If you have nerve damage
- If you have a skin disease or atopic dermatitis
- If you are a hospital outpatient and are taking medication
- If you have a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or if you have a concern about your skin when using the device

- **Do not use in the following cases or on the following areas:**

Failure to comply may cause accidents or other problems.

- If your skin is damaged, experiencing a rash, inflammation, or other abnormalities
- If you are under the influence of alcohol
- Head
- Anywhere inside the mouth, around the pubic region, or any mucous membrane
- Eyeballs
- Eyelids (eye sockets)
- Bulging moles
- Viral warts
- Areas where you have broken a bone
- Areas with internal bleeding, swelling, or scabbing following cosmetic medical procedures

*Avoid using the device for at least two weeks after undergoing any cosmetic medical procedure. After this period, consult a physician and only resume use if the physician confirms it is safe to do so (e.g., hyaluronic acid injections, silicone implants, Botox, thread lifts, double eyelid surgery, HIFU).

- Areas where metal was implanted by surgery (reconstructive surgery, plastic surgery, cosmetic surgery, etc.)

- **Do not use the device continuously or excessively, such as for commercial purposes.**

- **Do not place the device in water.**

Failure to comply may cause injuries or device malfunction.



Prohibition

WARNING (Use and storage)

- Do not use the device for any purpose or in areas other than those stated in the User's Manual.

Failure to comply may cause injuries or device malfunction.



Non-intended uses
are prohibited

- Keep the device and USB charging cable away from children and do not allow them to use either unsupervised. Do not allow children to play with the device.

- Children should obtain the consent of their guardians before using the device and should only use it under the supervision of a guardian or a specialist.

Failure to comply may cause accidents or injuries.

- Do not allow persons who are unable to express their intentions, physically challenged persons, or those whose lack of experience and knowledge prevents them from using the device safely on their own to use the device unattended.



User
restrictions

- Do not overuse the device.

For usage time and frequency, follow the instructions in the User's Manual.

Failure to comply may cause physical injuries.

- Before each use, check the device for any cracks, deformation, or other damage. Do not use if you see any sign of damage or failure.

Failure to comply may cause injuries or skin problems.



Prohibition

- Do not use, store, or charge the device in the following locations:

- In bathrooms or near bathtubs, showers, washbasins, or other containers filled with water
- Dusty areas
- Outdoors or a location where there is exposure to direct sunlight
- Uneven or unstable locations or a location where there is vibrations or shocks
- A location where there are strong radio waves or magnetic fields
- Near flammable substances (such as alcohol, thinner, benzene, nail polish remover, or aerosol sprays)
- A location where there is exposure to oil smoke or steam
- A location that may become hot, such as near a heating device (e.g., a heater)
- A location with drastic temperature changes (e.g., near an air conditioner or heater)

Failure to comply may cause fire, electric shock, accidents, device malfunction, falls, or injuries.

- Do not leave the device unattended while it is turned on.

Failure to comply may cause a fire or accidents.



Prohibition

WARNING (Use and storage)

- **Do not use if the USB charging cable is damaged.**

Failure to comply may cause electric shock, short circuit, or a fire.

- **Do not use the device in a dark room. In addition, do not look at, or directly into, the irradiation unit.**

Failure to comply may cause an eye injury.

- **Do not tighten with excessive force.**

Failure to comply may cause physical pressure, fatigue, and pain.

- **Do not use this product when asleep.**

Failure to comply may cause injuries or skin problems.



Prohibition

- **Do not disassemble, attempt to repair, or modify the device in any way.**

Failure to comply may cause a fire, electric shock, or other injuries.

Contact your local distributor for advice on repairing the device.



Disassembly prohibited

- **Do not dispose of the device in a fire or place it near fire or flames. Do not use, charge or store in a high-temperature or high-pressure area.**

This device contains a built-in rechargeable battery. Failure to comply may cause overheating, fire, or explosion.



Fire Hazard
Warning

- **When not using the product, or when cleaning it, make sure to turn off the power and to unplug it.**

When charging is complete, remove the charging USB charging cable from the main unit (controller section).

- **When unplugging the USB charging cable from the device and the USB port, do so by gripping the plug (do not touch the metal parts). Do not plug in or unplug the adapter with wet hands.**

Failure to comply may cause an electric shock, short circuit, fire, or device malfunction.

- **After use and when cleaning the device, be sure the power is switched off and disconnect the USB charging cable from the device.**

Failure to comply may cause fire or accidents.



Remove the
USB charging
cable

WARNING (Use and storage)

- If you feel an earthquake, see lightning or hear thunder, turn off the power and unplug the USB charging cable.

Failure to comply may cause fire or electric shock.



**Remove the
USB charging
cable**

- If the device is dropped or damaged, turn off the power immediately, pull out the USB charging cable, and request a repair/inspection.

- Do not use together with another electronic device (such as beauty equipment).

Failure to comply may cause skin problems or device malfunction.

- When using the device, remove all metal accessories (e.g., necklaces)

Failure to comply may cause skin problems.

- Ensure that the ambient temperature of the usage is properly maintained.

The range of ambient temperature in which the device can be used is from 10°C to 40°C. Verify the ambient temperature of the location where the device will be used. Failure to comply may cause the device to fail or malfunction.

The device may not correspond effectively to the ambient temperature of a new location immediately after moving to the new location.



Must Comply

- Do not apply excessive force to the mask.

Excessive force, including bending, pulling strongly, stepping on, or striking the mask, may cause damage or malfunction.



Must Comply

CAUTION (Use and storage)

- Do not press the device too hard and do not use it continuously.

Pressing the device stronger than necessary or using it continuously may damage the skin.

- Do not allow the device to come into contact with another beauty device or electronic appliance while the power is on.

Do not use with other electronic products (e.g., electric blankets/bed sheets, kotatsu tables) at the same time.

Failure to comply may lead to accidents due to malfunctions.

- Do not use while driving.

Failure to comply may cause accidents or injuries.



Prohibition

CAUTION (Use and storage)

- Use only Level 1 upon the first use and until you are accustomed to using the device. Even if you do not feel any stimulation, make only small adjustments to the intensity level. Use the device at an intensity level that suits you.

Make only gradual adjustments to the intensity level. Failure to comply may cause accidents or injuries.

- Always check for any signs of irritation on your skin during treatment.

If you experience any pain, fatigue, or skin abnormalities during use, stop using the device immediately.

If the device does not agree with your skin, stop using immediately. If symptoms persist the following day and beyond, consult a physician.

- If you are currently undergoing dental treatment or have teeth fillings, choose a lower intensity level and apply the device on a different area.

You may feel strong irritation or discomfort due to EMS on your cheeks or chin. If you feel strong discomfort, consult a physician.



Must Comply

- Avoid using the device on the same area for longer than 10 minutes.

When using the device on the same area again, be sure to leave an interval of 24 hours or more between treatments. Excessive use may damage the skin. The above is for reference only. Please observe your skin condition and use the device accordingly.



Must Comply

WARNING (The power supply)

- Use only with a 5 V DC USB power supply. Make sure the plug is fully inserted into the power outlet.

Failure to comply may cause a fire, electric shock, or device malfunction.

- If any problem occurs, turn off the power immediately.
- Always use the dedicated USB charging cable. Failure to comply may cause electric shock, short circuit, or a fire.



Must Comply

WARNING (The power supply)

- Do not use the device if the USB charging cable is damaged or if the USB port is loose.

Failure to comply may cause an electric shock, short circuit, fire, or malfunction. Contact your local distributor for an advice on repairing the device.



**Prohibition
of using a
damaged USB
charging cable**

- Do not connect too many plugs to one outlet exceeding the rated capacity of the outlet.

Failure to comply may cause an electric shock, short circuit, fire, or malfunction.



**Do not connect
too many plugs
to one outlet**

- Do not damage or modify the USB charging cable. Do not apply unnecessary force, bundle, or twist.

Failure to comply may cause the cable to break, resulting in a fire or electric shock. Straighten out any twists in the cord before storing.



**Improper
handling
prohibited**

WARNING (When using outside of Japan)

- Only charge the device using a 5 V DC USB power supply.

Failure to comply may cause fire or electric shock.



Must Comply

- Always check that the plug configuration and voltage match the specifications before use.

Failure to comply may cause overheating, fire, or malfunction.





**Fire Hazard
Warning**

- The plug configuration and voltage may vary depending on the region and the type of accommodation.
Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a consumer electronics retailer.


Names and Functions of Components

The design and specifications of the device and included accessories are subject to change for improvement purposes without notice.

Main unit (mask section)

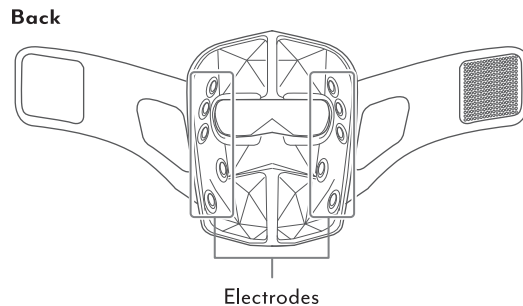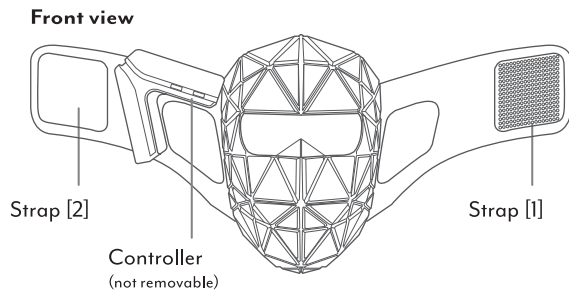



LED range

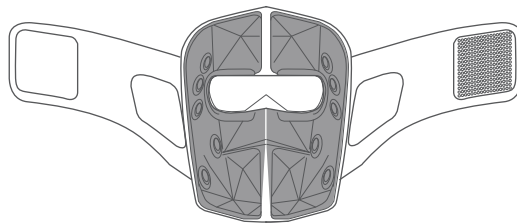

Main unit (controller section)

*Cannot be removed

EMS level / Charging / Remaining battery level indicators

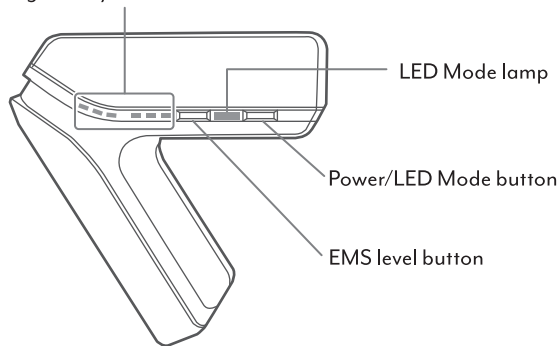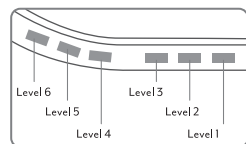

Accessories

● USB charging cable (Type-C)

● User's Manual

● Skincare product

Functions

Auto Power Off

The device will automatically switch off approx. 10 minutes after it has been turned on.


Mode/Level memory function

The device remembers the LED Mode and EMS level settings and returns to the same settings the next time it is switched on.


Before Use

Charge the Device

Using a household power outlet

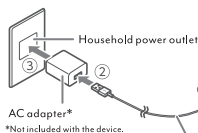

Using a computer

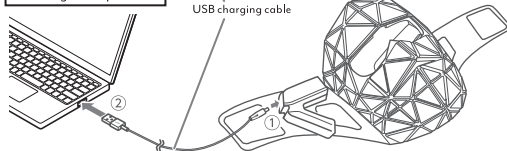

- Fully charge the battery before using the device for the first time after purchasing.
- When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be able to use the device for a short time only, even with a full charge. In such a case, charge the controllers again before using it.
- Use only the included USB charging cable.
- Please check the current rating indicated on the AC adapter (commercially available). Those with a rating of less than 2 A may not be able to recharge the device.

CAUTION

- When you have not used the device for a long period of time, check to see that it is functioning normally before restarting use. Failure to comply may cause accidents or malfunction.

Charging	The EMS level / Charging / Remaining battery level indicators for Levels 1–3 light up in sequence.
When charging is complete	The EMS level / Charging / Remaining battery level indicators for Levels 1–3 light up for 5 seconds then turn off.
Charging Time	Approx. 3 hours
Operating time	Approx. 50 minutes

*Actual charging time and operating time may vary depending on the remaining battery level and the charging conditions.

How to check battery level

You can check the remaining battery level by the blinking of Levels 1–3 of the EMS level / Charging / Remaining battery level indicators after powering on and the LED Mode lamp
→ EMS level / Charging / Remaining battery level indicators for Levels 1–6 light up in sequence.

Battery level	EMS level / Charging / Remaining battery level indicators		
	Level 3	Level 2	Level 1
High			
↕			
Low			
N/A*			

= Flashing = Lamp off



*When the charge runs out, a buzzer beeps six times, the LED Mode lamp flashes red three times, and the power turns off.

Pretesting on the skin

When using the device or the recommended skincare product for the first time, always perform a skin compatibility test to help prevent irritation or other skin problems.

- 1 Wash the inner side of your upper arm with a soap, then wipe dry with a towel**
- 2 Apply oil-free gel on washed skin and hold the controller on your upper arm. Make sure all three lower electrodes on either the left or right are in contact with your skin and select EMS Level 1 with your other hand, then use the device for approx. 3 minutes**

- For how to use, refer to the section “How to apply treatment” on page 16.

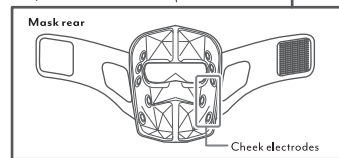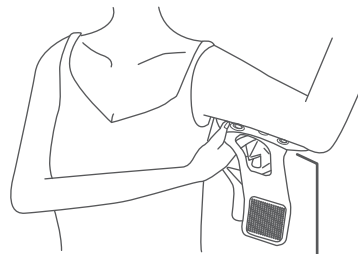

- 3 Wipe off the gel and check that your skin shows no irregularities**

- Use the device only if you experience no skin irregularities (such as redness or a burning sensation) within 24 hours.

CAUTION


- If you experience a rash, redness, itching, stinging, or other irritation at the tested part, stop the test immediately and wash the area thoroughly.
- Always check for any signs of irritation on your skin during treatment. If the device does not agree your skin, stop using immediately.

■ How to Use

Preparation

1 Wash your face and apply the included skincare product to the electrode parts of the device

*If you use the device on your skin when it is dry, it may be difficult for you to feel the EMS.

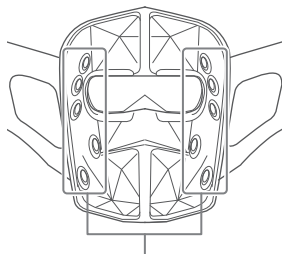

Electrodes



Care for your skin with skincare products.


- If you use your own skincare products, we recommend thick products.
- If your hair is long, we recommend tying it up.
- Cannot be used over a facial mask.

CAUTION

- Do not use oil-based skincare products, such as beauty oils or oil serums. It may cause deterioration of the mask.
- Refrain from applying flash treatments without first removing foundation and other makeup products that have been applied to the skin. Using the device without removing makeup may cause skin problems.

Wear the Blue Green Mask.

1 Fit the mask to your face

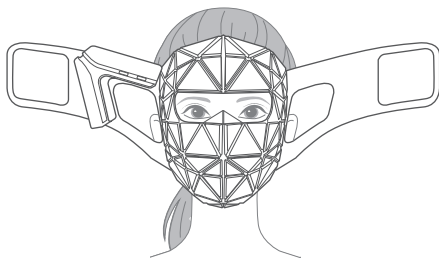

- First, align the chin line of the mask, then adjust and align it to the eye area.
- Make sure that no hair is caught in the mask.

2 Verify that the mask is firmly attached using a mirror



- Lay strap [2] on top of strap [1] and fasten while making adjustments so that the mask fits closely.
- Be careful not to fasten the straps too tightly.
- It is fine even if the mask fits closely to your mouth.

Tips


- If the mask is too large, tying your hair at the back of your head will help keep the face mask in place. If the mask is too small, you can use your hand to hold it in place without having to fasten the band.

<When using with EMS on>
Even if the mask does not fit perfectly and appears to be loose, as long as the section with all electrodes is in firm contact with your skin, you can proceed with the treatment.

<When using with EMS off>
Even if the mask does not fit perfectly and appears to be loose, as long as the section with the LED surface is in contact with your skin, you can proceed with the treatment.


How to apply treatment

1 Turn on the device

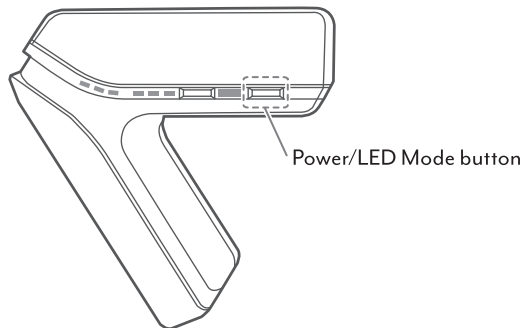

- Press and hold the Power/LED Mode button to turn on the power. After LED Mode lamp ⇒ EMS level / Charging / Remaining battery level indicators 1–6 light up in sequence, levels 1–3 of the EMS level / Charging / Remaining battery level indicators flash three times according to the remaining battery level (p. 14).

2 Choose the LED Mode (color)

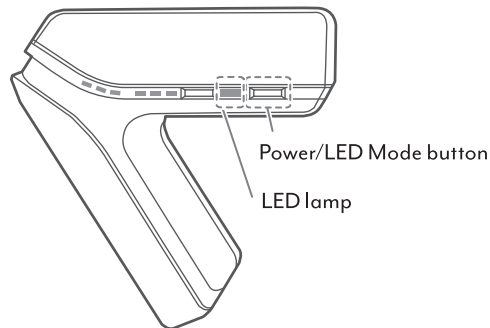


- Short press the Power/LED Mode button to switch the LED Mode (color).
- There are two colors for mode (color): GREEN and RED.
- The selected LED Mode (color) can be confirmed with the LED Mode lamp.

CAUTION

- Operate the controller while looking in a mirror.
- Refrain from using the device continuously for more than 10 minutes. Failure to comply may cause skin problems.
- If the light feels too bright during treatment, adjust the position of the mask or how the band is fastened or continue treatment with your eyes closed. Failure to comply may cause injury to your eyes or skin.

3 Select EMS intensity level

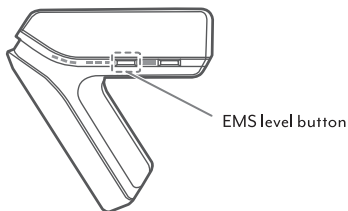


- Short press the EMS level button to switch between levels 0-6.
- Long press makes a “beep” buzzer sound, and the EMS level becomes 0 (Level 0 does not output EMS).
- You can check level switching with the EMS level / Charging / Remaining battery level indicators.

4 Turn off the device

- The device will automatically turn off approx. 10 minutes after being powered on.
- If you wish to stop mid-treatment, long-press the Power/LED Mode button on the controller.

*After each use, refer to “After Treatment” (page 20) and remove the remaining sebum, skincare product, and other particles that have adhered to the device.
Failing to remove them may cause device malfunction or skin problems.

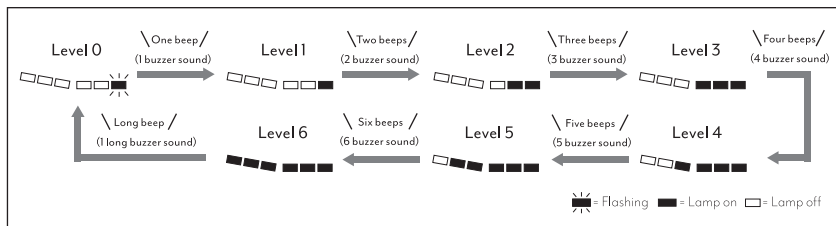

■ After Treatment

How to Clean

Clean the mask after each use.


1 Wipe off any dirt using a well-wrung soft cloth or similar materials

- If stubborn dirt remains, gently wipe with a cloth dampened with diluted neutral detergent.
- Do not leave the device with remaining sebum or skincare product on it.

2 Dry the device in a location that is well ventilated

CAUTION

- Always turn off the device before cleaning.
Failure to comply may cause electric shock, short circuit, or fire.
- Do not wash with water.
This device is not waterproof.
Failure to comply may cause device malfunction.
- Do not use alkaline cleaners, alcohol, thinner, benzene, or nail polish remover when wiping.

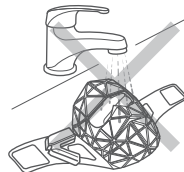

Storing the Device

- Store the device in a clean, well-ventilated place out of the reach of children.
- Remove the USB charging cable from the controller when not in use for a long period of time.
- Do not leave the device in hot or humid places, such as near open flames, in direct sunlight, or in a bathroom.

Maintenance and Inspections

- If you have not used the device for an extended period, inspect it carefully and check that the device works properly and safely before use.
- If the device is malfunctioning, go to 'Troubleshooting' (page 21~22) and try to identify the problem before arranging to have the device repaired. If the problem persists, contact your local distributor.
- If the device malfunctions or appears damaged, please contact your local distributor.

Additional Information

Troubleshooting

If you notice an abnormality during use, check for the following symptoms.

Symptom	Possible causes and solutions
The device does not turn on	● Press and hold the Power/LED Mode button for about two seconds (page 18).
	● The battery may be depleted. Charge the battery (page 14).
The device only operates for a few minutes even after charging	● The device will automatically turn off approx. 10 minutes after being powered on.
	● The battery level is low. Charge the battery (page 14).
	● The battery may have reached the end of its life.
	● The battery may not be charging properly. Check that charging has started by seeing whether the EMS level / charging / remaining battery indicators are lit or flashing.

■ Additional Information

Troubleshooting

If you notice an abnormality during use, check for the following symptoms.

Symptom	Possible causes and solutions
The LED Mode lamp flashes during operation	<ul style="list-style-type: none">● The battery level is low. Charge the battery (page 14).
The device does not start charging	<ul style="list-style-type: none">● Make sure the AC adapter (commercially available) and USB charging cable are properly connected (page 14).● Check the current rating indicated on the AC adapter (commercially available). Adapters with a rating lower than 2 A may not be able to charge the device.
EMS feels too strong (tingling or painful)	<ul style="list-style-type: none">● The EMS level may be set too high. Please lower the level. Always start at level 1 when using the device for the first time.
EMS feels weak (no sensation)	<ul style="list-style-type: none">● The EMS level may be set too low. Please increase the level.
No changes in EMS intensity are felt / Tingling or pain is felt even after adjusting the EMS level	<ul style="list-style-type: none">● Your skin may be too dry. Apply lotion, gel, water, or similar products to moisturize the skin.
EMS sensation varies by treatment area	<ul style="list-style-type: none">● Wipe off any dirt on the electrodes.
	<ul style="list-style-type: none">● The electrodes are not properly in contact with the skin. Ensure the mask is worn correctly.
	<ul style="list-style-type: none">● The battery level is low. Charge the device (page 14).
	<ul style="list-style-type: none">● Depending on muscle mass, some areas may be less sensitive to EMS.
An electronic sound is heard during treatment	<ul style="list-style-type: none">● This is a normal operating sound and does not indicate a malfunction.

Frequently Asked Questions

Q

What is the recommended duration and frequency of use?

A

The device can be used every day. However, as a general guide, use it up to 10 minutes per session, once per day.

*The device will automatically switch off approx. 10 minutes after it has been turned on.

*The above information is a guideline only. Please check your skin condition while using the device.

*Always check for any signs of irritation on your skin during treatment.

If the device does not agree your skin, stop using immediately.

Q

Can I use the device in the bathroom?

A

Do not use the device in the bathroom. This device is not waterproof.

Q

Can I use the device while sleeping?

A

Do not use the device while sleeping.

Q

Can I use the device on dry skin?

A

It can be used when the EMS level is set to 0.

When using EMS on dry skin, it may be difficult to feel the EMS stimulation.

If the EMS feels weak, please moisten your skin using the included skincare product before use.

Q

Can I use the device over a facial mask?

A

Not recommended. The LED and EMS effects do not work properly.

Q

My skin has developed a rash, inflammation, or itchiness

A

This may be caused by stimulation from heat, light, or EMS. If the symptoms persist the following day or beyond, consult a physician.

Q I am concerned about the smell of the mask.

A A characteristic odor from the silicone rubber may be noticeable, but this is not abnormal.

Q I feel some discomfort or fatigue.

A This may be caused by stimulation from heat, light, or EMS. Stop using the device until the discomfort in your face subsides.

Q My teeth hurt when I use the device.

A The pressure of the mask or EMS may cause tooth pain. If you feel discomfort, consult a physician.

Q Can I use the device immediately after suntanning?

A Using the device immediately after suntanning may cause inflammation. If any skin problems occur, stop using the device immediately and consult a physician.

Q My skin has developed redness.

A The pressure from the mask may cause temporary redness. If the redness persists the following day, consult a physician.

Q The mask leaves a mark on my face.

A The mask is made of stretchable silicone. Depending on your skin condition, it may leave temporary marks. However, there is no need for concern, as the marks will fade over time.

Q Is there a possibility that makeup color will transfer to the mask?

A Yes, color transfer may occur. Please remove makeup before using the device.

Product Specifications

Rated Power Supply	DC5V 2A
Power Consumption	Maximum approx. 8.0 W (when charging)
Battery	Lithium-ion battery
Charging Time	Approx. 3 hours (Actual charging time may vary with conditions of the charging environment)
Operating time	Approx. 50 minutes (from a fully charged state) *Operating time may vary with conditions of use.
Outer Dimensions	With mask open: W635 × D112 × H210 (mm) With mask folded: W348 × D60 × H245 (mm)
Weight	Approx. 520 g
Materials	Mask: silicone and stainless steel Controller: ABS resin
Country of Origin	Made in China
Accessories	One USB charging cable (Type-C) User's Manual One skincare product
Recommended ambient temperature	10°C–40°C

■ After-sales Service

Recommended skincare products

*Sales may be discontinued without prior notice.

Flawless Serum Hydrating



QTY: 80 g

Item No.: CS21004

Rechargeable battery

- The rechargeable battery is a consumable item (It is not covered by warranty).

The battery is expected to last 300 charge cycles. If the device no longer operates for more than a short period of time even after the battery has been recharged, the battery may have reached the end of its service life.

For battery replacement, please contact customer service.

Disposing of the device

- Follow your local government rules for proper disposal of this product.
- The battery must be removed from the device before it is disposed of.
- The lithium-ion battery used in this product is a valuable recyclable resource.



Removing the Battery


- Be sure to turn the power off before removing the battery.
- Use a screwdriver to remove the battery, being careful not to injure yourself.
- When removing the battery, be careful not to damage it with a hard object. Failure to comply may cause a generation of heat or a fire.
- Follow your local government rules for proper disposal or recycling of this product.

Operation Steps


Refer to figures ① to ④ below to remove the battery.

- ① Remove the four screws on the controller cover.


*If they are difficult to remove, expose the screw parts by lifting the strap.

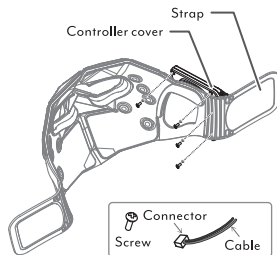

- ② Release the four latches on the controller cover and open the cover.

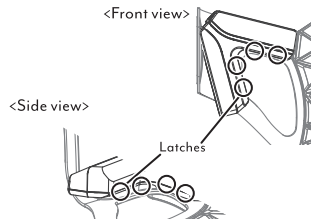


- ③ Disconnect the three cable connectors.

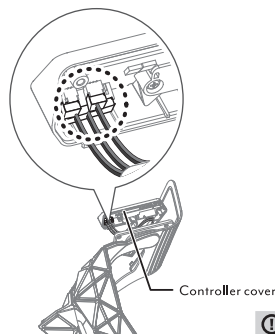

- ④ Disconnect the connector of the cable attached to the battery, insert a thin object like a plastic ruler under the battery, and remove the battery.

*It is attached with strongly adhesive tape.

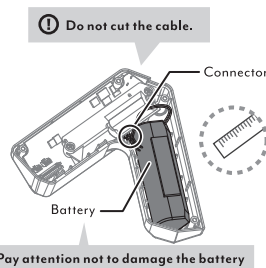

The figures above are intended for product disposal purposes only and not for repair.

CAUTION: Disassembling or repairing the device yourself may cause malfunction.

Contacts for customers outside Japan

- Please contact your local distributor or visit YA-MAN Global website (<https://www.ya-man.com/en/>).

Marking Description



This marking indicates that this product should not be disposed with other household wastes throughout the UK and EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased, they can take this product for environmental safe recycling.



This marking indicates that it is a detachable supply unit.

This product must be charged only using the included detachable power supply.

Memo

Memo

