

User's Manual

Household Beauty Device Medi Lift Eye EPE-10

- Read the 'Safety Precautions' on the other side of this sheet carefully to use the device safely.

©2019 YAMAN LTD. ALL RIGHTS RESERVED.

Unauthorized reproduction or reprinting is prohibited.

EPE10-2-001E

Introduction

- Thank you for purchasing this YAMAN product.
- Read this User's Manual carefully to use the device properly and safely.
 - After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
 - Make sure to include this User's Manual if you are transferring ownership of the device to another person.

Accessories

- AC adapter
- Magnetic USB charging cable*
- User's Manual


*Do not use these accessories with any products other than Medi Lift.

The design and specifications of the device and peripheral accessories are subject to change for improvement purposes without notice.


Names and Functions of Components

Controller

●Front

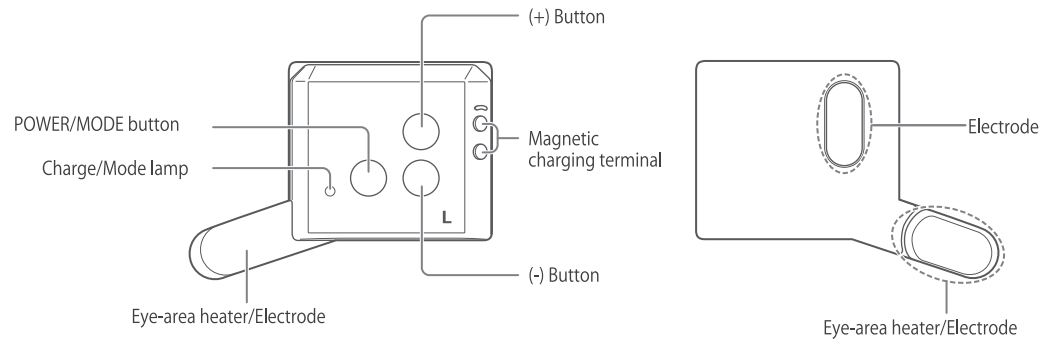


●Rear

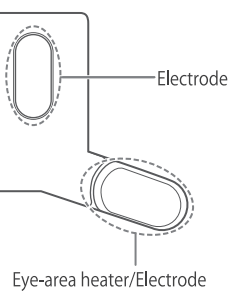


Silicone mask

●Front

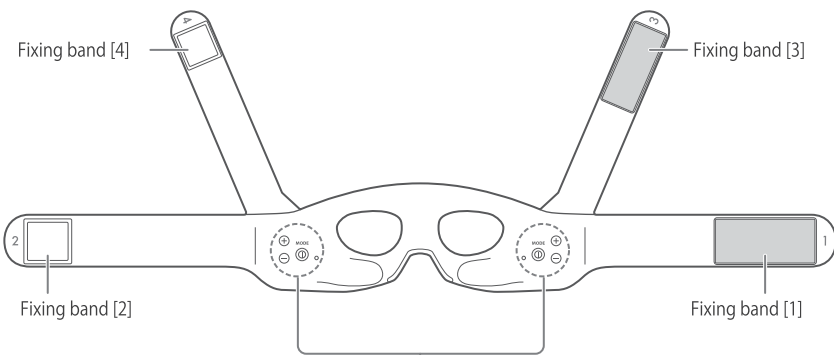

Controller operating parts

●Rear

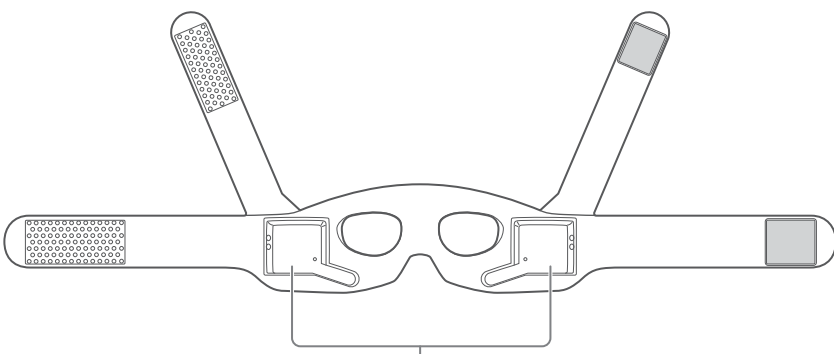

Controller attachment parts

Functions

Auto Power OFF

The device will automatically switch OFF approx. 10 minutes after it has been turned ON and the mode was set.

If the mode is changed during use, the device will switch OFF approx. 10 minutes after this change, so please be mindful of how many minutes the device has been in use.

EMS

Stimulates muscles electrically.

Heater function

The heater function warms the eye area.

Micro-current

The weak electrical current conditions the skin.



Before use

Charge the device

How to remove the controller

Hold onto the eye-area heater part of the controller and lift upwards to remove it.

- Be careful not to use excessive force when removing the controller. This may damage the silicone mask or main unit electrode parts.

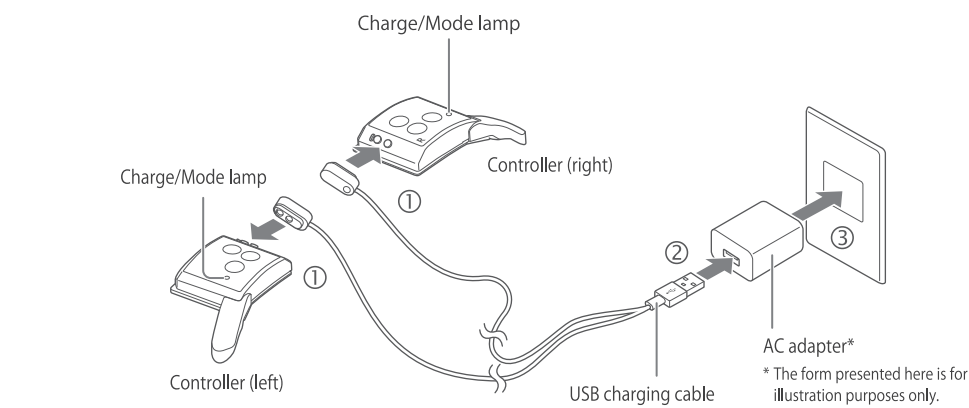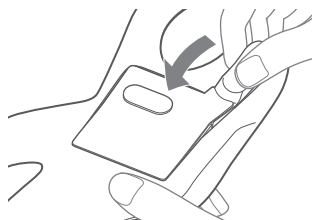

- Charging will begin when the controller is connected as shown in the diagram. While charging, the Charge/Mode lamp will blink in red.
- When charging is complete, the Charge/Mode lamp will stop blinking and remain red, turning OFF after one minute. Remove the USB charging cable.
- At an indoor temperature of approximately 25°C, it will take about 2.5 hours to reach a full charge, which will allow the device to operate for about 3 hours. Actual charging time and operating time will differ depending upon the remaining power and the charging environment.

ATTENTION

- Be sure to attach the magnetic charging terminal in the correct direction. Charging will not begin if the magnetic charging terminal is attached in the wrong direction.
- If there is moisture on the controller or the magnetic charging terminal, be sure to wipe the water off thoroughly. If there is moisture on the magnetic charging terminal while charging, the device may malfunction.
- When charging the battery for the first time after purchasing the device, use the device only after the battery is fully charged.
- When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge displayed. In such a case, charge the device again before using it.
- Use only the enclosed USB charging cable and AC adapter.

Pre-testing on Skin

When using the device or a skin-care product for the first time, always test for compatibility with your skin.

ATTENTION

- Be sure to perform a pretest in order to prevent skin problems from occurring.
- If a rash, reddening, itchiness, stinging sensation, or other abnormality occurs at the pre-testing area, immediately stop the pre-test and wash the pre-testing site.
- Do not use the device if any skin abnormality is perceived.

- 1 Wash the inner side of upper arm with soap and water, then wipe dry with a towel.
- 2 Apply your regular eye-area serum to the washed area of your upper arm.
- 3 Hold the device against your upper arm with the opposite hand so that two electrodes are touching the skin, then operate the device on Lift Mode Level 1 for approx. 3 minutes.
 - Test the device with the controllers attached to the silicone mask. Hold the device so that the two electrodes on the back of the silicone mask on either the left or right side are touching the skin.
 - When the power is turned ON, the device is set to Lift Mode and the EMS level (intensity) is set to 0 (MIN).
 - Refer to 'TREATMENT' for information on how to use the device.
- 4 Wipe off the eye-area serum and check that there is no skin abnormality.
Continue to use the device only if you experience no abnormalities (reddening or burning sensation) on your skin within 24 hours.

How to Use

Check out the instruction video

An instruction video of demonstrating how to use the device is available. Scan the QR code on the right side.

<https://bit.ly/2YuRNXO>



Wear the silicone mask on your face.

Attach controllers (2) to the silicone mask.

Be mindful of the direction of the controllers when attaching them.

- 1 **Rear**
Left Right
Check the direction of the eye-area heater part of each controller.
- 2 **Align the position of one controller with the grooves in the silicone mask and press the controller gently into the mask.**
- 3 **Press the controller down and attach it firmly to the silicone mask. Attach the remaining controller to the silicone mask in the same way.**
- 4 **Check that the eye-area heater part of each controller is slightly raised.**

ATTENTION

- Be careful not to use excessive force when attaching each controller to the silicone mask. Doing so may cause damage to the controller.
- Do not press the controller buttons when attaching each controller to the silicone mask.
- Be sure that each controller is turned OFF while you are attaching it to the silicone mask.

Place the silicone mask on your face.

After washing your face to remove makeup, please thoroughly moisturize your skin using your regular skin-care routine. We recommend using Medi Lift Eye Serum for moisturizing the eye area.

The EMS will not conduct an electric current if there is no moisture between your skin and the electrodes on the silicone mask. Be careful to ensure that there is always moisture between your skin and the electrodes.

ATTENTION

- Please do not use beauty oil or any other cosmetics containing oil. Doing so may cause the silicone mask to deteriorate.
- If you put the silicone mask on your face without removing foundation, the silicone mask and/or the controllers may be stained. Be sure to remove all makeup before putting on the silicone mask.

- 1 Tie your hair back and pin up your bangs to ensure that your forehead is not caught between your skin and the silicone mask.
- 2 Turn the silicone mask over and apply an appropriate amount of skin-care product to the four electrodes (two on each side) on the controllers.
 - Be sure to apply the skin-care product AFTER you have attached the controllers to the silicone mask.
- 3 **Align the silicone mask with the position of your eyes.**
 - Position the mask and make sure that the eye-area heaters on the controllers are not in contact with the eyeballs.
 - Align the nose position of the mask with your nose and place the mask on your face, ensuring the word 'YAMAN' on the front of the mask is not upside-down.
- 4 Pull fixing bands [1] and [2] diagonally upwards and attach them firmly with the fastener.
 - Head shapes differ from person to person. Fix the fixing bands in positions where they are held firmly in place and will not slip.
 - If the fixing bands slip down, take measures such as tying your hair in a ponytail to keep the fixing bands firmly in place.
- 5 Pull fixing bands [3] and [4] upwards to the top of the head and attach them firmly with the fastener.
 - Adjust the fitting of fixing bands [3] and [4] by adjusting their attachment angle.
 - The silicone mask can also be used without fixing bands [3] and [4].
- 6 Look in a mirror to check that the silicone mask is firmly attached to your face.
 - Depending on facial shape, the silicone mask may not fit properly or there may be space between your skin. However, it can be used as long as the four electrodes (two on each side) are firmly in contact with the skin.

Cannot feel EMS?

Is there any difference in EMS sensation between the left and right sides of the face?

- The skin lacks moisture.
Thoroughly moisturize the skin around your eyes using your regular skin-care routine. We recommend using Medi Lift Eye Serum.
- When you wish to increase the EMS stimulation, apply an ample amount of basic skin-care product directly to the four electrodes on the controllers. Differences in moisture levels between the left and right sides of the face may cause a difference in sensation. EMS sensation between the left and right sides of the face may also differ depending on the degree of contact between the electrodes and the skin.


ATTENTION

- Be sure to put on the silicone mask AFTER you have attached the controllers.
- Be sure that the controllers are turned OFF while you are putting on the silicone mask.
- Be sure to tie your hair back when putting on the silicone mask.
- Be careful not to get your hair caught in the fasteners on the silicone mask when you are putting on the mask.
- Check to make sure that the four electrodes on the back of the silicone mask (two on each side) are firmly in contact with the skin. If the contact points are loose, it may not be possible to use the device correctly.
- If you use an oil-based leave-in hair conditioner, wait until your hair becomes dry before using the device.
- Be careful to ensure that there is always moisture, such as eye-area serum, between your skin and the electrodes. Insufficient moisture may not only prevent the device from functioning correctly, but also reduce sensation and cause skin problems.

Treatment

Until you become familiar with how to use the silicone mask, please look in a mirror while you are using it.


- 1 **Turn the power ON**
Press the POWER / MODE buttons on the left- and right-hand sides of the controller simultaneously for at least one second.
A very long beep will sound followed by the Charge/Mode lamp light in purple.
 - Be sure to turn the power ON when you are wearing the silicone mask with the controllers attached.
 - When the power is turned ON, the device is set to Lift Mode and the EMS level (intensity) is set to 0 (MIN).

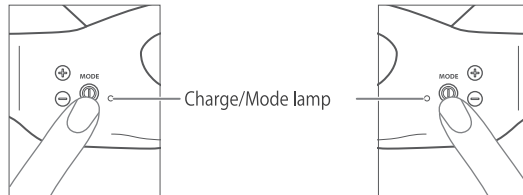

The device will automatically switch OFF approx. 10 minutes after the mode was set.

If you want to stop the treatment during use, press both left and right POWER / MODE buttons simultaneously for at least one second to turn the power OFF.

- 2 **Select your preferred mode.**
Press the POWER / MODE buttons on the left and right controllers simultaneously. The buzzer will sound each time you press the button and switch the mode.

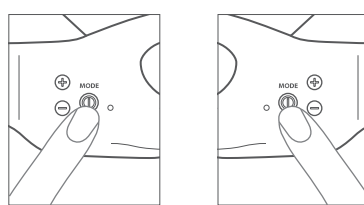


Lift mode (purple)
Provides 'tapping' stimulation.

→

Stretch mode (blue)
Provides 'tightening' stimulation.

- 3 **Select your preferred EMS level (intensity).**


Press the +/- buttons on the left and right hand controllers simultaneously.

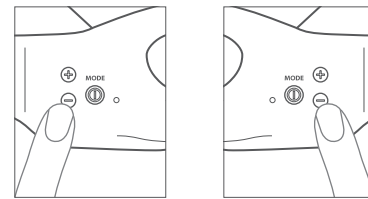

- The EMS level (intensity) of each mode can be adjusted to five different levels.
- On Level 0, only the heater can be used.
- The heater level cannot be changed. The heater takes approx. one minute to heat up.
- As shown below, a different sound will beep for each setting.

Level 0	One long beep
Level 1	One beep
Level 2	Two beeps
Level 3	Three beeps
Level 4	Four beeps
When pressing the minus (-) button with the device on Level 0*	One very long beep
When pressing the plus (+) button with the device on Level 4*	One very long beep

* When pressing the minus (-) button with the device on Level 0 or when pressing the plus (+) button when the device is on Level 4, one very long beep will sound, indicating the MAX/MIN intensity.

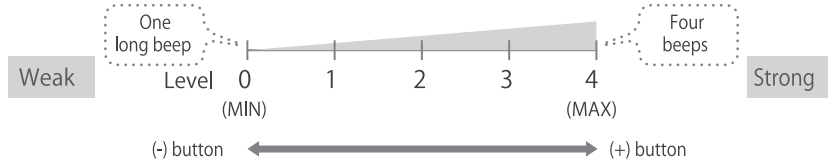

ATTENTION

- Use only Level 1 upon first use and until you become familiar with how to use the silicone mask.
- If your skin becomes dry, please apply additional eye-area serum during the treatment. Please turn the device OFF and remove the silicone mask before applying the additional serum to the eye area.
- If you should feel an unpleasant tingling pain, please do not increase the EMS level (intensity).

END

The device will automatically switch OFF approx. 10 minutes after the mode was set. When the power switches OFF, the device emits two long beeps.

Rub any remaining eye-area serum into your skin.

If you want to stop the treatment during use, press both left and right POWER / MODE buttons simultaneously for at least one second to turn the power OFF.

As a general guide, use the device once a day (Maximum of twice a day).

- The device will automatically switch OFF approx. 10 minutes after the mode was set.
- The above information is a guideline only. Please observe your skin's condition and use the device accordingly.
- If you perceive any skin abnormality, immediately stop using the device.
- When using the device continuously, we recommend having intervals of five hours or longer between treatments.


After Treatment

How To Clean (Controllers)

To keep the product clean, please be sure to clean the controllers after each use. Remove the controllers from the silicone mask.

WARNING

- Do not clean the device, be sure to turn the power OFF.
May cause an electric shock, short circuit, or fire.
- Do not wash in water.
This device is not waterproof.
May cause a failure of the device.

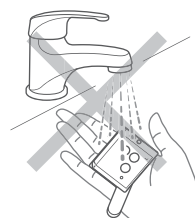

- 1 **Wet a piece of soft cloth with water and firmly wring it dry before wiping away stains.**
For stubborn stains, wipe gently with a cloth soaked in a diluted neutral detergent.

ATTENTION

- Do not leave gel or eye-area serum on the device.
- Refrain from using an alkali cleaning agent, alcohol, thinner, benzene, nail polish remover, or other such substance.
- Due to the characteristics of the device, strongly rubbing the electrode parts may cause color staining. When maintaining the device, be sure to use a cloth that you do not mind getting dirty.




How To Clean (Silicone mask)

To keep the product clean, please be sure to clean the silicone mask after each use.

Remove the controllers from the silicone mask.

- 1 **Remove the controllers and wash the silicone mask in warm water.**
For stubborn stains, wipe gently with a cloth soaked in a diluted neutral detergent.
- 2 **Use a towel or dry cloth to wipe off moisture.**
Be sure to wipe the device thoroughly to ensure that no moisture or eye-area serum remains on the electrode.
- 3 **Dry the device in a location with good ventilation.**
Leaving the device with grime or moisture on it could cause mold, discoloration, or odors.

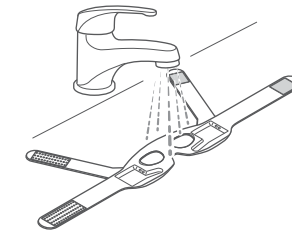

ATTENTION

- Do not leave eye-area serum on the device.
- Do not use hot water.
- Refrain from using an alkali cleaning agent, alcohol, thinner, benzene, nail polish remover, or other such substance.



Safety Precautions

- Read these Safety Precautions carefully before using the device, and use it correctly.
- Sensibly use the device in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any malfunction or unusual reaction to treatment occurs.
- Do not use the device in a defective state, as doing so could cause physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or ATTENTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

●Classifications of Precautions

DANGER

Improper use of this device poses a severe and likely risk of death or serious injury.

WARNING

Improper use of this device may pose a risk of death or serious injury.

ATTENTION

Improper use of this device may result in injury or property damage*.

* 'Property damage' refers to damage affecting a dwelling, household belongings, or domestic animals and pets.

●Symbols

This symbol shows you what you must not do.

This symbol shows you what you must do.

This symbol shows you what you must be careful with.

DANGER (regarding use)

- The following people should not use the device.
Use of the device could cause an electronic medical device to malfunction.
 - ◆ Do not use in conjunction with any medical electronic equipment.
 - pacemaker or other implanted electronic medical devices;
 - artificial heart-lung machine or other electronic life-support medical devices, and
 - electrocardiograph or other worn electronic medical devices.

WARNING (regarding use and storage)

- Do not use under the following conditions:
May cause an accident or other problems.
 - If you are under a doctor's care •While growing or during adolescence •Acute disease
 - Tuberculous disease •Fever •Hypertension •Infectious disease •Malignant tumor •Blood disease
 - Visceral disease •Skin disease •Allergic constitution •Atopic dermatitis •If you have skin problems.
 - Sensitive skin •If you have unhealthy skin blemishes caused by physical irritation
 - If you have been diagnosed with loss of temperature sensation
 - If you have an illness of the eyes (cataracts, glaucoma, retinal detachment, fundal hemorrhage, myodesopsia, etc.)
 - Excessive myopia (±60 or higher) •If you had an eye or face operation
 - If you have diseases of the eyes or face or you are under treatment
 - If you are severely fatigued •If your physical condition is unwell
 - If you are under the influence of alcohol •While pregnant or nursing or during menstruation
 - If you take medication that induces drowsiness
 - If you have capillary dilatation due to long-term use of steroid hormonal drugs or a liver disorder
- Do not use the device in the various areas indicated below or outside of specified areas on the face.
May cause an accident or other problems.
 - Cuts •Areas where surgery (including plastic surgery) was performed
 - Anywhere inside the mouth or any mucous membrane •Irritated or red skin caused by sunburn •Acne
 - Skin with purulent inflammation •Places where skin is itchy or feels hot •Bulging moles •Viral warts
 - Adam's apple •Eyes •Eyelids

* If you are uncertain as to whether it is safe to use the device, consult a physician prior to use.
- Do not use the device continuously for an overly long period (for example, through commercial use).
May cause an injury or a failure of the device.
- Do not immerse the controllers in water.
May cause an injury or a failure of the device.
- Do not use the device near a flammable substance (such as alcohol, thinner, benzene, nail polish remover, or spray).
Device may explode or catch fire.
- Do not leave the device unattended while it is turned ON.
May lead to the device catching fire or other similar accidents.

- Do not use the device on infants or people who are unable to express their will. Use by unattended people with physical disabilities should be avoided. A minor should use the device only under the supervision of his or her guardian who has given their consent.
May cause an accident or injury.
- Keep away from children and do not allow them to use the device. Do not allow children to play with the device.
May cause an accident or injury.

- Do not use the device for any purpose other than those stated in the User's Manual.
May cause an injury or a failure of the device.

- Do not use the device when an odor or abnormal generation of heat is detected.
Failure to comply may cause damage or a fire.

- Do not disassemble, attempt to repair, or modify this device in any way.
Contact your local distributor for advice on repairing the device.

- Do not place the device in an unstable location. Do not drop, bump, forcefully push, or knock the device, or poke it with a sharp object.
May cause an injury, accident, failure of the device, electric shock, or fire.

- Do not sit or stand on the device, or place heavy objects on it.
May cause an injury or a failure of the device.

- During use, if your eyes or around the eye areas feel hot, stop using the device immediately. Do not attempt to tolerate the heat.
Doing so may lead to low-temperature burns.

ATTENTION (Use and storage)

- Do not use anywhere on the body except for the face.
May cause an injury.
- Do not use the device on pets.
May cause an injury.
- During treatment, carefully check that your skin is not showing any abnormality.
Do not use the device if the treatment does not agree with your skin.
- Before use, check to make sure that the electrode is not cracked, deformed, or broken.
Do not use the device if it is damaged or malfunctioning.
May cause your skin to suffer an injury or other problems, or lead to existing conditions worsening.
- Do not wipe the device with alcohol, thinner, benzene, nail polish remover, or other such substance.
Failure to comply may cause damage to the surface of the device.
- Do not leave the device water or in the bathroom or any other high-humidity environment.
Do not use the device in a sauna or spa.
May cause an injury or a failure of the device.
- Never recharge the device in the bathroom or near a receptacle containing water.
May cause an electric shock, short circuit, fire, or failure of the device.
- Do not use together with another electronic device (beauty equipment).
Failure to comply may lead to skin troubles or malfunctions.

* This product can be used together with some of other Medi Lift series. Read User's Manual of your purchased carefully product before using the device.

- Do not use the device in outdoors.
May cause a failure of the device.

- When using the device continuously, we recommend having intervals of 5 hours or longer between treatments.
May cause an injury or generation of heat.
- If you should begin to feel ill or a skin abnormality appears during use, stop using the device immediately.
If symptoms continue the following day and beyond, consult a physician.
- Use only Level 1 for the EMS level (intensity) upon first use and until you become familiar with how to use the silicone mask. Even if you do not feel the EMS stimulation, make only small adjustments to the intensity levele.
May cause an accident or injury. Be sure to check the condition of your skin while using the device.
- If you feel an unpleasant tingling on your skin, reduce the EMS intensity level.
Failure to comply may cause skin damage.
- After each use, refer to 'After Treatment' and remove eye-area serum and other particles that have adhered to the device.
May cause a failure of the device or skin problems.
- Be sure to comply with the following items.
May cause an accident, injury or a failure of the device.
 - Avoid using the device for longer than the set time.
 - Use the device appropriately based on your skin conditions.
 - Do not use the device if your skin is experiencing a rash, inflammation, itchiness, or other abnormalities.
 - Refrain from using the device together with your skin routine immediately after sun-tanning/being sunburnt or peeling.
- Be careful to ensure that there is always moisture, such as eye-area serum, between your skin and the electrodes.
If your skin is dry, reapply the eye-serum after taking the silicone mask off with the power turned OFF.
If the level of moisture is low, EMS flow may decline, reducing effects and causing skin problems.

- Make sure to remove makeup before using the device.
If you use the device over makeup, the foundation may cause dark staining of the skin.
- During treatment, carefully check that your skin is not showing any abnormalities.
Stop using the device immediately if the treatment does not agree with your skin.
- When you have not used the device for a long time, check it to see that it is functioning normally before restarting use.
May lead to accidents or malfunctions.
- When using the device, remove all metal accessories (watch, necklace, rings, earrings, etc.)
Failure to comply may cause skin problems.
- Stop using the device and turn it off if you feel a tremor or hear thunder.
Disconnect the AC adapter when the device is being charged.
- While using the device, if you feel that a contact lens has shifted, stop use immediately and remove the contact lens.

WARNING (The power supply)

- Do not damage, modify, bundle, or otherwise apply undue force to the AC adapter or USB charging cable.
Do not place heavy items on top of the AC adapter or USB charging cable, or pinch them between other items. May cause damage, fire or electric shock.
- Untwist cords before storage.
May cause a disconnection or electric shock.
- Do not move or lift the controller by pulling the AC adapter or USB charging cable.
May cause an electric shock, short circuit, fire, or disconnection.

- Do not use if the AC adapter or USB charging cable is damaged or if the wall socket is loose.
May cause an electric shock, short circuit, or fire.
Contact your local distributor for advice on repairing the device.

- Do not wet, wash, or immerse the AC adapter, USB recharging cable, or controllers in water.
May cause an electric shock, short circuit, fire, or failure of the device.

- Be sure to turn off the power when not in use or when cleaning the device. Remove the AC adapter and USB charging cable from the controller and the power outlet when not in use. Do not insert or disconnect the plug with wet or damp hands.
May cause an electric shock, short circuit, or fire.

- When disconnecting the AC adapter and USB charging cable from an outlet, be sure to pull it out by holding the plug instead of the cord (Do not touch the metal prongs). Do not insert or disconnect the plug with wet or damp hands.
May cause an electric shock, short circuit, or fire.

- Use with an operating voltage range from 100 V to 240 V.
Make sure that the AC adapter and USB charging cable are fully inserted into an electrical outlet. Be sure to check the plug configuration and the voltage of the power outlet before using the device in countries other than Japan.
May cause a fire, electric shock, or a failure of the device.

- In the unlikely event that an abnormality arises, turn the device off immediately.
Refrain from using the device if faulty operation, smoke, abnormal noises, or odors are detected, or if it is otherwise operating in an abnormal state.

- Be sure to use the included AC adapter and USB charging cable with this product. Plug the AC adapter directly into a wall outlet and do not overload an electric circuit.
May cause an electric shock, short circuit, fire, or failure of the device.
 - * The USB charging cable included with the device is for recharging use only.

- Before charging the device, check to make sure that there is no serum or water droplets on the charging terminal.
May cause an electric shock, short circuit, or fire.

- Be sure to turn off the power when not in use or when cleaning the device.
Make sure to disconnect the AC adapter after completion of charging.

- Do not leave the device near fire, in the hot sun, in a bathroom or other high-humidity or high-temperature location, or near combustible gas or flammable objects. Do not dispose of the device in a fire or place it near fire or flame. Do not use or charge in a high-temperature area or if the device is hot. Do not bring the device close to a heater.
The device features a built-in rechargeable battery. May cause the generation of heat, a fire, the explosion of the device, or failure of the device.

Additional Information

Storing

- Store in a clean, well-ventilated place out of the reach of children.
- Remove the AC adapter and USB charging cable from the controller and the power outlet when not in use for a long period of time.
- Do not leave the device in hot or humid places, such as near flame, in bright sunlight, or in a bathroom.

Maintenance and Inspections

- If you have not used the device for an extended period, inspect it carefully and check that the device works properly and safely before using it again.
- If the device is malfunctioning, go to 'Troubleshooting' and attempt to identify the problem before arranging to have the device repaired. If the problem persists, contact your local distributor with a proof of purchase.
- If the device malfunctions, appears damaged, or needs service, please contact your local distributor.

Troubleshooting

If you perceive any abnormality while using the device, verify the following points before arranging to have the device repaired.	
Problem	Possible causes and solutions
Does not turn on.	→ Press the Power/Mode button for about one second.
	→ Insufficient charge. Charge the battery.
The device only operates for a few minutes even after charging.	→ The device will automatically switch OFF approx. 10 minutes after the mode was set.
	→ The battery is low. Charge the battery.
	→ The battery has expired. Contact local a distributor.
The device does not start charging.	→ Connect the AC adapter to an outlet.
	→ Check to make sure that the magnetic charging terminal is connected in the correct direction.
EMS is too strong (EMS causes unpleasant tingling or pain)	→ Lower the EMS intensity with the (-) button.
	→ Your skin is dry. Moisturize your skin by applying the eye-area serum.
EMS is weak (EMS cannot be felt).	→ Increase the EMS intensity with the (+) button.
	→ The electrodes are not placed on the skin correctly. The silicone mask is not attached correctly.
The heater does not warm up (Warm-up is very slow.)	→ Your skin is dry. Moisturize your skin by applying the eye-area serum.
	→ It takes approx. one minute for the heater to warm up
	→ The battery is low. Charge the battery.

Frequently Asked Questions

Q What is the recommended duration and frequency of use?	A As a general guide, use the device up to 10 minutes at one time. * The device will automatically switch OFF approx. 10 minutes after the mode was set. * The above information is a guideline only. Please observe your skin's condition and use the device accordingly. * If you perceive any skin abnormalities, immediately stop using the device.
Q Can I use the device constantly?	A We recommend using this device once a day, for 10 minutes. When using the device continuously, we recommend having intervals of five hours or longer between treatments.
Q How much eye-area serum should I use each time?	A This device is designed to achieve better EMS effects with sufficient use of eye-area serum. You may feel a slight tingling sensation.
Q Can I use the device in the bathroom?	A This device is not waterproof. Do not use the device in the bathroom.
Q Can I use the device while sleeping?	A Do not use the device while sleeping.
Q There is a difference in stimulation between the right and left sides of my face.	A Because the faces are asymmetrical, you may feel different stimulation on the left and right sides of your face. Adjust the right and left EMS levels.
Q I can hear a buzzing noise in my ear.	A That is the sound of the electricity flowing, not a failure of the device.
Q My skin is experiencing a rash/inflammation/itchiness.	A This is caused by the EMS stimulation. If symptoms continue the following day and beyond, consult a physician.
Q I feel some discomfort and tiredness.	A This is caused by the EMS stimulation. Do not use the device for a while until the feeling of discomfort on your face goes away.
Q My teeth hurt when I use the device.	A If your teeth are being treated or you have a silver crown, your teeth may hurt when you are using the device. If you feel discomfort, consult a physician.
Q Can I use the device while pregnant?	A Do not use the device during pregnancy or immediately after childbirth. Consult your family physician.
Q Can I use the device immediately after suntanning or shaving?	A Using the device after suntanning or shaving may cause inflammation. If skin problems should occur while using the device, stop use immediately and consult a physician.
Q My skin is experiencing redness.	A The heater or pressure from the silicone mask may cause temporary redness. If the redness continues the following day, consult a physician.
Q The mask leaves a mark on my face.	A The mask is made of elastic silicon and its purpose is to physically lift the skin using the elastic band. Therefore, it is designed to place pressure on the face. Depending on the condition of your skin, the mask may leave a mark. However, there is no need for concern as the mark will disappear with time.
Q The belt is too short and tightens my skin.	A Extension bands are available for purchase. If you wish to purchase an extension band, contact the local distributor.

Product Specifications

Rated power supply	DC5V 0.5A
Power consumption	Maximum approx. 3.0 W (during recharging)
Battery	Lithium ion battery
Charging time	Approx. 2.5 hours (May vary depending on the charging environment.)
Operating time	Approx. 3 hours (from fully charged state) * Operating time may vary with conditions of use and EMS level.
Outside dimensions	Silicone mask: approx. W 620 × H 230 (mm) (including attached bands when lying flat) Controller: Approx. W78 × H54 × D15 (mm) (including electrodes)
Weight	Silicone mask: approx. 98 g Controller: Approx. 32 g each
Materials	Silicone mask: silicone Controller: ABS resin
Country of origin	Silicone mask: Made in China Controller: Made in Japan
Accessories	AC adapter..... One Magnetic USB charging cable One User's Manual
Recommended environmental temperature	10 °C to 40 °C

After-sales Service

Rechargeable battery

- The rechargeable battery is a consumable item. (It is not covered by the product warranty).
- If the device no longer operates for more than a short period of time even after the battery has been recharged, the battery may have reached the end of its service life (The battery is expected to last 300 charge cycles).

Disposing of the device

- Follow your local government rules for proper disposal of this product.
- This product uses a lithium ion battery.



Consultation counter for warranty from outside Japan

For countries and regions outside Japan, please contact the store where you purchased for warranty details.
Or check through YAMAN Global Official Website <https://www.ya-man.com/en/>