



User's Manual

Home-care Beauty Device

Photo Plus Deep Lift

YJFA1

Features


This product uses RF radio frequency function (RF) to warm the skin. It is a facial beauty device with ionic cleansing, ion moisturizing, EMS, Photo, micro-current, and cooling functions.

Users can choose a combination of various modes to carry out care procedures according to the needs of individual skin conditions.

- RF...The heat reaches deep into the skin, so that the skin can be cared for in a more suitable state
- Ion derivation...Suctions dirt which cannot be removed with a regular cleansing with ion derivation
- Ionphoresis...Unique wave form of iontophoresis helps easier penetration of lotion
- EMS...Able to stimulate muscles with micro-current
- LED...Skin care with LED light irradiation
- Micro-current...Care for around the eyes with a weak electric current
- Cool down...Apply to the skin surface with ice-cold temperature to make the skin firmer
- Vibration...Care with fine vibrations

※to stratum corneum

Names of Components

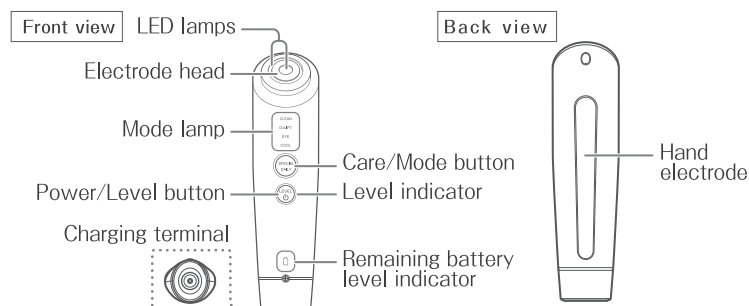

Accessories



- **Charging base**



- **Cotton stopper**
(attached to the electrode head of the main body when packaged)
Use this stopper to fix a cotton to the electrode head

- **USB charging cable (USB Type-C)**

- **User's Manual**

*Exclusively designed for this product. Do not use attachments other than those provided with the product.

*The design, specifications of the device and peripheral accessories are subject to change for improvement purposes without notice.

Modes

	Mode	Function	Cotton pad	Recommended cosmetics* ¹	Operating time	Level	Recommended frequency
Daily Care	DAILY	<ul style="list-style-type: none"> •RF •EMS •Iontophoresis •LED (Red) 	Either with or without a cotton pad	Skin lotion or serum	Approx. 3 minutes	6 levels	Every day* ²
Special Care	CLEAN	<ul style="list-style-type: none"> •RF •Ion conduction •Vibration (only during ion emission) 	With a cotton pad	Skin lotion	Approx. 3 minutes	6 levels	2 to 3 times a week
	D×LIFT	<ul style="list-style-type: none"> •RF •EMS •DYHP •LED (Red) 	Either with or without a cotton pad	Serum or sheet mask	Approx. 5 minutes	6 levels	
	EYE CARE	<ul style="list-style-type: none"> •RF •EMS •Micro-current 	Without a cotton pad	Serum for around the eyes	Approx. 2 minutes	6 levels	2 to 3 times a week* ³
	COOL	<ul style="list-style-type: none"> •Cool down •LED (Blue) 	Without a cotton pad	Place directly on the skin* ⁴	Approx. 3 minutes	-	Every day* ⁵

*1: If using your regular skincare product, be sure it is water soluble. Do not use oil-based or scrubbing skincare products. Failure to comply may cause product malfunctions or skin problems.

*2: Do not use on the same day as Weekly Care Mode (except in EYE CARE Mode and COOL Mode).

*3: If you are concerned about dryness of the skin around the eyes, we recommend using the device 5 times a week.

*4: We recommend using the device after moisturizing your skin with a milky lotion or other skincare products.

*5: The above recommendations are merely guidelines. Please adjust the frequency of use, such as by using the device 2 to 3 times a week while observing the condition of your skin.

Read the instructions below carefully before use for a proper usage.

Precautions

- Be sure to perform a pre-test before using the device for the first time.
- Be sure to move the device on your skin while it is in use. Using the device continuously on the same area may cause burns or other injuries. (except in EYE CARE Mode and COOL Mode).
- The warmth of the RF (radio frequency) changes depending on the type of skincare products used together with the skin's condition at the time of treatment.
- Do not use the device over the upper eyelids (except in COOL Mode), lower eyelids or outer corner of the eye (except in EYE CARE Mode and COOL Mode). *The LED light can be irradiated to the upper eyelids, lower eyelids and outer corner of the eye. Close your eyes and hold the device at least 5 cm away from your skin. Do not look directly into the LED light.
- When holding the device, make sure that you are holding the hand electrode.

Basic Method of Use

- 1 Remove makeup and gently cleanse your face**
Remove oil and dirt from your skin so that it is in a clean condition prior to a treatment.
- 2 Attach cotton pad to the electrode head if necessary and apply skin lotion to the cotton pad**
- 3 COOL Mode can be used over make-up. Turn on the device**
Hold down the Power/Level button (hold down the button for approx. 3 seconds or longer). The device will emit a beeping sound, the Care/Mode button (DAILY) lights up, and DAILY Mode starts.
- 4 Select the mode**
To change the mode, press and hold or briefly press the Care/Mode button. If Daily Care is selected, the Care/Mode button (DAILY) lights up. If SPECIAL Care is selected, the mode button (SPECIAL) and the mode lamp for the selected mode light up.
- 5 Select the treatment level (intensity) (except in COOL Mode)**
Each time you briefly press the Power/Level button, the device will beep once and the level will change.
- 6 Hold the device with your hand touching the hand electrode (the silver rim) on the rear**
If your hand is not touching the hand electrode, it cannot operate correctly.
- 7 Apply the electrode head to your skin and move the head over the area being treated**
 - When DAILY Mode is finished, the device will automatically turn off.
 - When the treatment in each SPECIAL Care Mode is finished, the device will go into Standby (except when in COOL Mode). Select your preferred mode to start the next mode.
 - If you want to turn the power off during use, press and hold down the Power/Level button.
- 8 Care for the skin after treatment**
When you have finished treatment, rub the remaining lotion into your skin before using your usual basic skincare products to condition it.

How to Use the Device on Your Skin


- Recommended speed for moving the device: approx. 2 cm per second
- If the treatment temperature feels too hot, move the head more quickly.
- One short beep will sound to indicate the half-way point in the operating time.
When you hear the beep, begin care on the other side of your face.



Areas you should not treat


DAILY Mode

Move the electrode head over the skin in the direction of the arrows shown on the right. The LED light can be used to treat the upper eyelids, lower eyelids and outer corner of the eye.

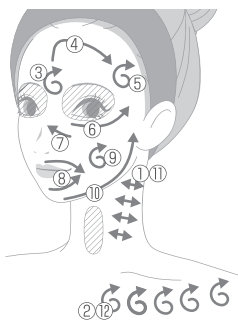


CLEAN Mode

Move the electrode head over the skin in the direction of the arrows shown on the right.

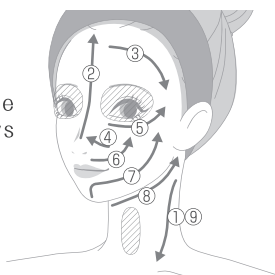

D×LIFT Mode


Move the electrode head over the skin in the direction of the arrows shown on the right. The LED light can be used to treat the upper eyelids, lower eyelids and outer corner of the eye.

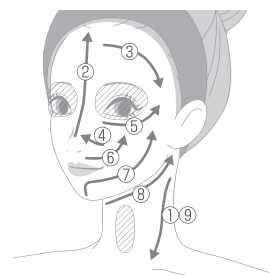

EYE CARE Mode

Press the electrode head repeatedly on your face from the inner corner to the outer corner of the eye as if stamping on the skin (Press → Release → Move electrode head). Press down on the skin for approx. 3 seconds. If you are concerned about dryness of the skin around the eyes, we recommend using the device for 2 minutes per treatment, twice a day, 5 times a week.

- * Do not press down continuously on the skin for longer than 3 seconds.
- * Do not move the device while pressing.

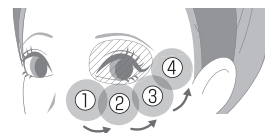

COOL Mode

Press the electrode head repeatedly on your face as if stamping on the skin (Press → Release → Move electrode head).

- * This mode can also be used over the upper eyelids, lower eyelids, and outer corner of the eye.
- * When using the LED light on the eyelids, be sure to close your eyes and keep the electrode head at least 5 cm away from the skin.