


How to Use

Household Optical Beauty Device Rei Beauté Cool Power YJEA9N


Names and Functions of Components

Main Unit

*This manual explains using illustrations when the roller attachment is attached.
*The xenon lamp is a consumable item.

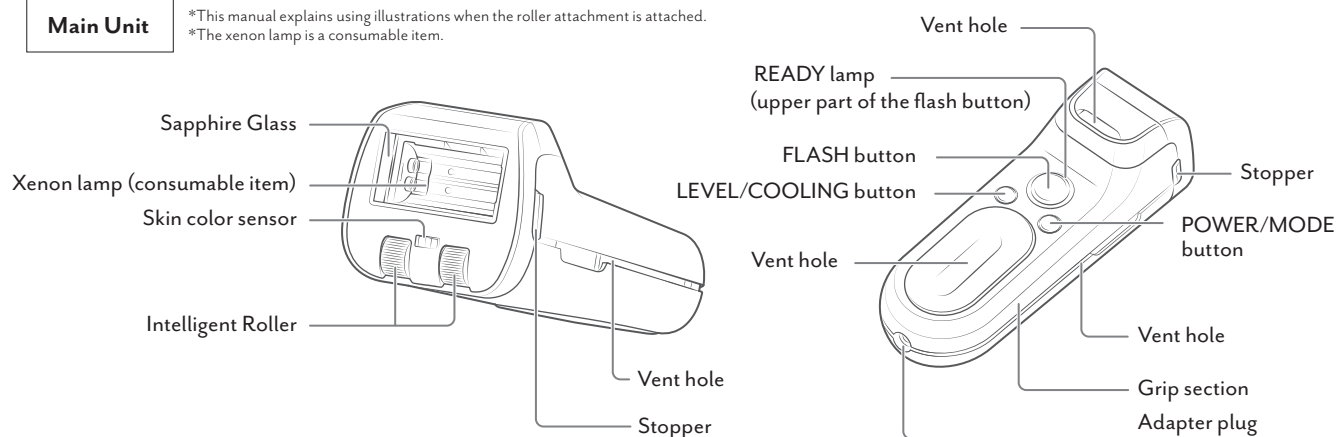

Accessories

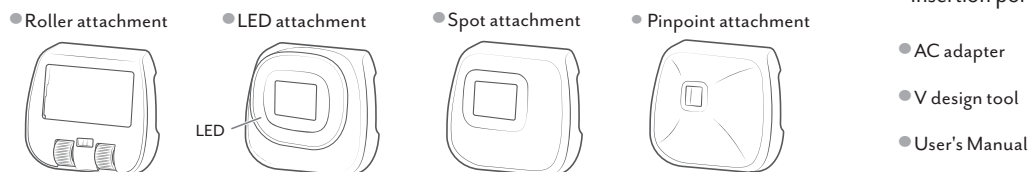

Modes

Select the mode and level

You can change the treatment mode by pressing POWER/MODE button.

Press the LEVEL/COOLING button to select a treatment level. Mode and level will switch as below.

Roller attachment

[Mode]
Single Mode → Roller Mode

[Level]
1 (MIN) → 2 → 3 → 4 → 5 (MAX) → Auto Level (AUTO)

[Cooling feature]
Cooling ON → Cooling OFF

*Press and hold the LEVEL/COOLING button to change ON and OFF.

LED attachment


[Mode]
Single Mode → Blue-green LED Mode

[Level] • Single Mode
1 (MIN) → 2 → 3 → 4 → 5 (MAX)

[Level] • Blue-green LED Mode
1 (MIN) → 3 → 5 (MAX)
Continuous Flash Single Flash

Spot attachment / Pinpoint attachment

[Mode] Single Mode
[Level]
1 (MIN) → 2 → 3 → 4 → 5 (MAX)

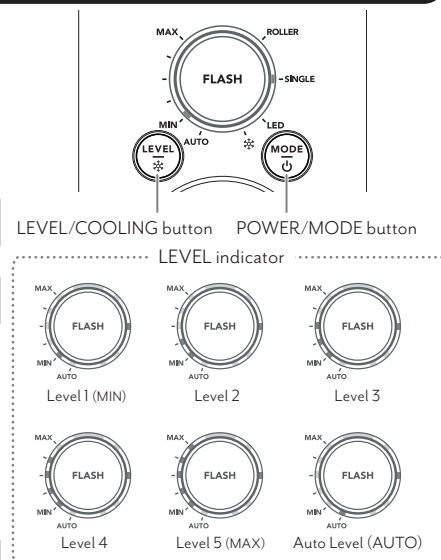


Please refer to the table below for available modes, levels, and functions for each attachment.

		Roller attachment	LED attachment	Spot attachment	Pinpoint attachment
Mode	Roller Mode	○	—	—	—
	Single Mode	○	○ (Red LED)	○	○
	Blue-green LED Mode	—	○ (Blue-green LED)	—	—
Level	Auto Level (AUTO)	○	—	—	—
Function	Cooling feature	○	—	—	—

Single Mode

(Apply to all attachments)


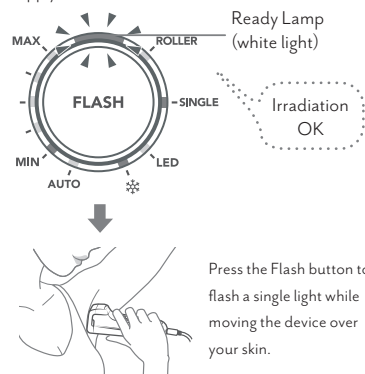


Recommended areas for use


Roller attachment	LED attachment	Spot attachment	Pinpoint attachment
- Underarms - Bikini line - Face (cheeks)	- Face (cheeks) - Face (around the mouth) *The device cannot be used on the lips. *When using the flashing light on the skin below the nose, stretch the skin to prevent the light from flashing on the lips. *Close your eyes when you bring the LED close to your face.	- Fingers - V-line - I-line - O-line *Do not flash the light inside the I-line, such as on male genitalia.	- Face (brow area) - Face (forehead) - Face (nose) *Do not use on eyelids or around the eyes. *If the light feels too bright during treatment, close your eyes. *When using on the brow area, wear sunglasses or protective goggles to protect your eyes.

Roller Mode

(For a Roller attachment only)

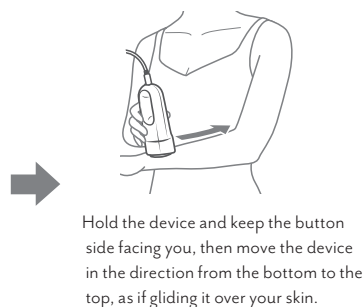
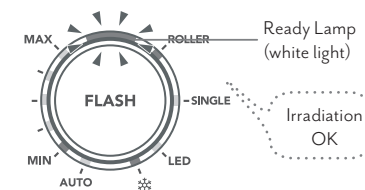


Recommended areas for use

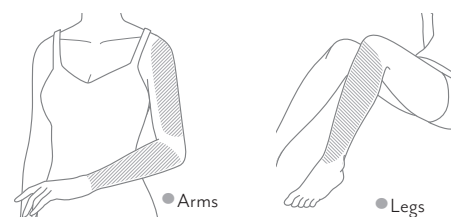

Blue-green LED Mode


(Only available for LED attachment)

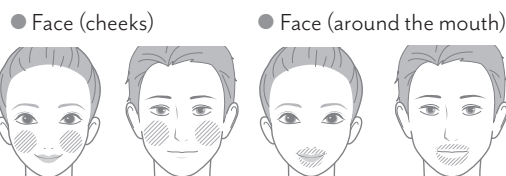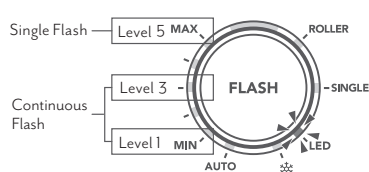

During continuous flashing (up to 5 flashes), glide the device gently over your skin.
For single flashing, press the flash button and move the device, flashing one area at a time.
* If the LED light feels too bright, close your eyes during use.

Cooling feature

(For a Roller attachment only)

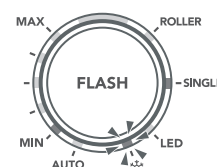

Press and hold the Level/Cooling button to change the sapphire glass cooling function ON and OFF.

How to use

1

Prepare the treatment area


Remove unwanted hairs from the area to be treated with an electric shaver or razor.



2

Attach the attachment and turn on the power


Hold down the POWER/MODE button.

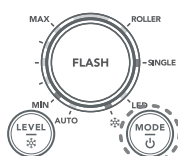

3

Select the mode


Press the POWER/MODE button.

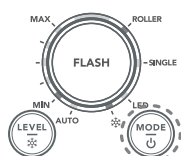

4

Select the level

Press the LEVEL/COOLING button.

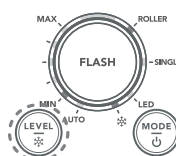

5

Flashing light

Press the device firmly to your skin and start a treatment.



ATTENTION

- Be sure to use the attachment that is appropriate to the area being treated. Failure to comply may cause a burn injury.
- Ensure that neither you nor anyone else looks directly into the treatment window of the device. Also ensure that no one nearby looks at the treatment window. Failure to comply may damage your eyes.
- Moisturize your skin daily to prevent your skin from becoming dry. If your skin is dry, the device might not work properly.
- When treating the face, press the treatment window firmly against the skin. If the treatment window is not pressed firmly against the skin, the flashing light may become overly bright.