



How to Use

Household Optical Beauty Device Rei Beauté Cool Pro YJEA6W

Names and Functions of Components

Main Unit

*This manual explains using illustrations when the roller attachment is attached.
*The xenon lamp is a consumable item.

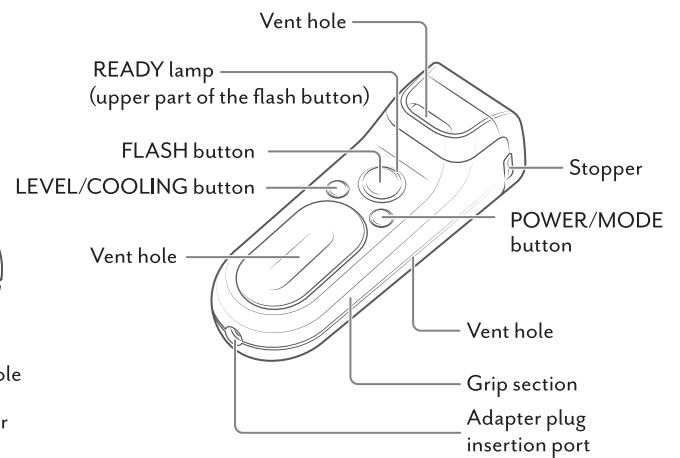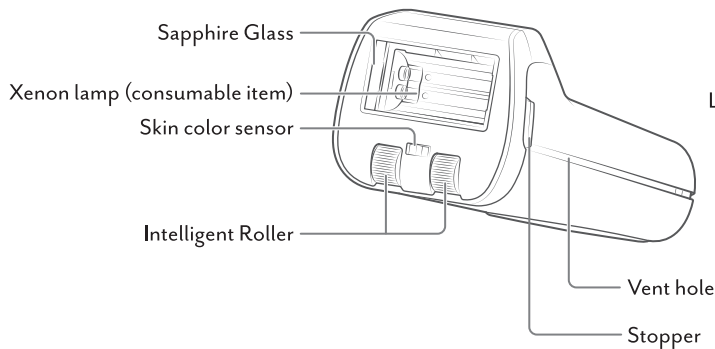


Accessories


● Roller attachment



● LED attachment

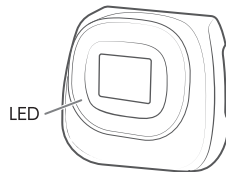

● Spot attachment

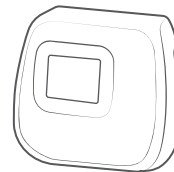

● AC adapter

● V design tool

● User's Manual

Modes

Select the mode and level

You can change the treatment mode by pressing POWER/MODE button.

Press the LEVEL/COOLING button to select a treatment level. Mode and level will switch as below.

Roller attachment

[Mode]

→ Single Mode → Roller Mode*

*Roller Mode and Auto Level are not available when the LED attachment or the spot attachment is mounted.

[Level]

→ 1 (MIN) → 2 → 3 → 4 → 5 (MAX) → Auto Level (AUTO)

[Cooling feature]

→ Cooling ON → Cooling OFF

*Press and hold the LEVEL/COOLING button to change ON and OFF.

LED attachment

[Mode]

→ Single Mode → LED Mode

[Level]

→ 1 (MIN) → 2 → 3 → 4 → 5 (MAX)

*There is no treatment level selection available for LED Mode.


Spot attachment

[Mode]


Single Mode

[Level]

→ 1 (MIN) → 2 → 3 → 4 → 5 (MAX)

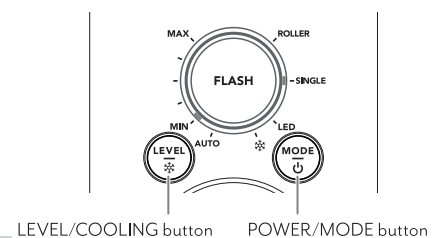


LEVEL/COOLING button POWER/MODE button

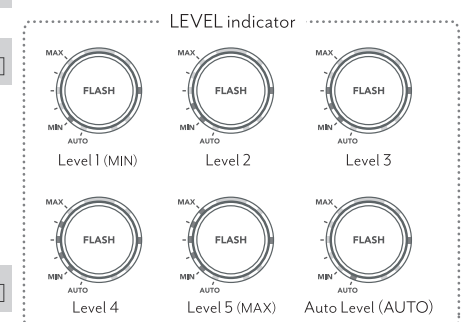

Single Mode

(Apply to all attachments)

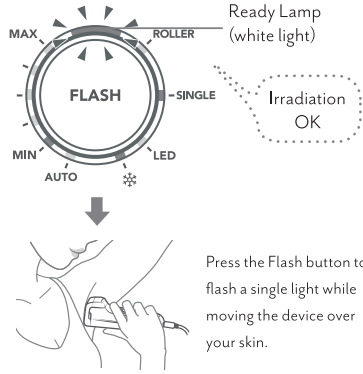


Recommended areas for use


Roller attachment	LED attachment	Spot attachment
<ul style="list-style-type: none"> ● Underarms ● Bikini line 	<ul style="list-style-type: none"> ● Face (cheeks) ● Face (around the mouth) <p>*The device cannot be used on the lips. *When using the flashing light on the skin below the nose, stretch the skin to prevent the light from flashing on the lips. *Close your eyes when you bring the LED close to your face.</p>	<ul style="list-style-type: none"> ● Fingers ● V-line ● I-line ● O-line <p>*Do not flash the light inside the I-line, such as on male genitalia.</p>

Roller Mode

(For a Roller attachment only)

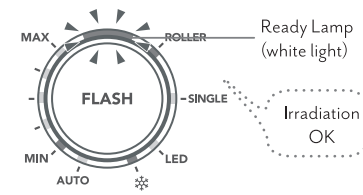


Recommended areas for use

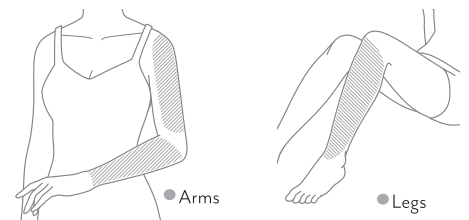


LED Mode

(Only available for LED attachment)

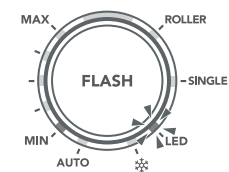


Press the treatment window to your skin as if stamping on it.

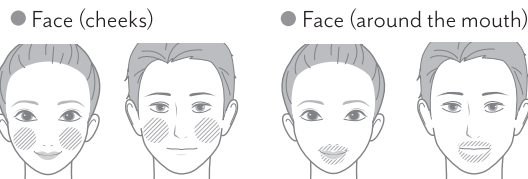

Recommended areas for use

Cooling feature

(For a Roller attachment only)

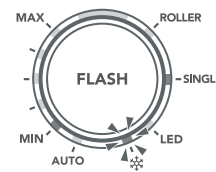

Press and hold the Level/Cooling button to change the sapphire glass cooling function ON and OFF.

How to use

- ### 1 Prepare the treatment area

Remove unwanted hairs from the area to be treated with an electric shaver or razor.
- ### 2 Attach the attachment and turn on the power

Hold down the POWER/MODE button.
- ### 3 Select the mode

Press the POWER/MODE button.
- ### 4 Select the level

Press the LEVEL/COOLING button.
- ### 5 Flashing light

Press the device firmly to your skin and start a treatment.

- Be sure to use the attachment that is appropriate to the area being treated. Failure to comply may cause a burn injury.
- Ensure that neither you nor anyone else looks directly into the treatment window of the device. Also ensure that no one nearby looks at the treatment window. Failure to comply may damage your eyes.
- Moisturize your skin daily to prevent your skin from becoming dry. If your skin is dry, the device might not work properly.
- When treating the face, press the treatment window firmly against the skin. If the treatment window is not pressed firmly against the skin, the flashing light may become overly bright.

ATTENTION