

User's Manual

Household Optical Cosmetic Device | Rei Beautè Go | STA-207

Table of Contents

Introduction.....	1	After Treatment.....	25
Safety Precautions	1	Additional Information	27
When Using Outside of Japan	11	● Troubleshooting.....	27
Names and Functions of		● Error Indications.....	29
Components	12	● Frequently Asked Questions	30
Modes	15		
Before Use	17	● Product Specifications	31
How to Use	19	After-sales Service.....	32



Introduction

Thank you for purchasing this YA-MAN product.

- Read this User's Manual carefully in order to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Be sure to include this User's Manual if you are transferring ownership of the device to another person.

Safety Precautions

- Read these Safety Precautions carefully before using the device in order to use it correctly.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any malfunction or unusual reaction to treatment occurs.
- Do not use the device in a defective state, as doing so could cause physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, and ATTENTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

●Classifications of Precautions



DANGER

Improper use of this device poses a severe and likely risk of death or serious injury.



WARNING

Improper use of this device may pose a risk of death or serious injury.



ATTENTION

Improper use of this device may result in injury or property damage*.

* 'Property damage' refers to expanded damage affecting a dwelling, household belongings, or domestic animals and pets.

●Symbols



This symbol shows you what you must not do.



This symbol shows you what you must do.



This symbol shows you what you must be careful with.

WARNING (Power supply)

- **Always use the dedicated AC adapter.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.



- **Do not use the device if the AC adapter is damaged or if the plug is loose.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.
Contact the store where you have purchased for advice on repairing the device.



- **Do not move or lift the device by pulling the AC adapter.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.



- **Do not wet, wash, or immerse the AC adapter in water.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.



- **In the unlikely event that an abnormality arises, immediately unplug the AC adapter.**

Refrain from using the device if faulty operation, smoke, abnormal noise or odors are detected, or if it is otherwise operating in an abnormal state.



- **When disconnecting the AC adapter from an outlet, pull it out by grasping the plug instead of the cord (Do not touch the metal prongs). Do not insert or disconnect the plug with wet hands.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.

- **When the device is not in use or when cleaning the device, make sure to disconnect the AC adapter from the outlet.**

Failure to comply may cause an electric shock, short circuit, or fire.




- **Do not overload an electric circuit.**

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.



- **Do not damage, modify, bundle, or otherwise apply undue force to the AC adapter.**

Placing a heavy object on or pinching the electrical power cord could damage the cord or cause a fire or electric shock.

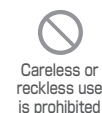

- **Untwist cords before storage.**

Placing a load on the cord may cause disconnection or electric shock.

- **Do not wrap the cord around the device or AC adapter when storing the device.**

Placing a load on the power cord, which may result in disconnection, electric shock, or short circuit.

- **Only use at AC100-240V in order to avoid a potential fire or electrical shock hazard. Make sure that the AC adapter is fully inserted into a power outlet.**

Failure to comply may cause a fire, electric shock or accident.



- **When not using the device for an extended time, disconnect the AC adapter from the outlet.**

Failure to comply may result in a fire.

- **After use, be sure to turn off the power and disconnect the AC adapter from the outlet.**

Failure to comply may cause the device to catch fire or may lead to an accident.

DANGER (Use and storage)

● The following people should not use the device.

Use of the device could cause an electronic medical device to malfunction.

◆ Do not use in conjunction with any medical electronic equipment.

- Pacemakers or other implanted electronic medical devices; · Artificial heart-lung machines or other electronic life-support medical devices, and
- Electrocardiographs or other worn electronic medical devices.



WARNING (Use and storage)

● Do not use under the following conditions:

Failure to comply may cause an accident or other problems.

- If you are under a doctor's care · During adolescence · If you have an acute disease · If you have tuberculous · If you have a fever
- If you are suffering from hypertension · If you have an infectious disease · If you have a malignant tumor · If you have a febrile disease
- If you have a blood disease · If you have a visceral disease · If you have a heart disease · If you have a skin disease
- If you have a melanoderma · If you have an allergic constitution · If you have atopic dermatitis · If you have hypersensitivity to light
- If you have skin problems · If you have sensitive skin · If you have unhealthy skin blemishes caused by physical irritation
- If you have skin irritation or redness caused by sunburn · If you are suffering from dermatitis due to cosmetics · If you are severely fatigued
- If your physical condition is unwell · If you are under the influence of alcohol · While pregnant or nursing or during menstruation
- If you take medication that induces drowsiness · If you have capillary dilatation due to long-term use of steroid hormonal drugs or a liver disorder;
- If you had an eye or face operation; · If you have diseases of the eyes or face or you are under treatment; · Excessive myopia (-6D or higher).



● Do not use the device on the various areas indicated below or outside of specified areas on the face.

Doing so may cause an accident or other problems.

- Cuts · Areas where surgery (including plastic surgery) was performed
- Anywhere inside the mouth, lips, around the pubic region, or any mucous membranes · Acne · Areas of the skin with purulent inflammation
- Areas with itchiness or hot sensations · Bulging moles · Viral warts · Adam's apple · Over implants · Eyeballs · Eyelids
- Areas around the eyes · Sunburnt skin · Head · Bone of the nose · Right above the heart (around the stomach) · Epilepsy
- Injuries (scars or other dark-colored areas) · Blemishes, moles · Eyebrows, the forehead, hair on the head, the nape of the neck · Ears
- Nipples, areolae · Inside the navel · Areas with makeup · Areas where a large vein or other blood vessels can be seen
- Tattoos or body-painted areas

*If you are uncertain as to whether it is safe to use the device, consult a physician prior to use.

WARNING (Use and storage)

- **Do not use in the vicinity of children.**

As the flash constitutes a strong source of light, injury to the eye may result.

- **Do not use the device when the vent holes are blocked.**

Failure to comply may result in smoke, fire, or burns.

- **Do not touch the filter (glass plate) immediately after using the flash.**

As the filter becomes very hot immediately after use, touching it with your finger may cause a burn injury.

- **Do not use the device for overly long periods (such as for commercial purposes).** Ten minutes of use should be followed by at least ten minutes of rest.

Failure to comply may cause an injury, the generation of heat, or the development of cracks in the glass tube.

- **Do not use the device near a flammable substance (such as alcohol, thinner, benzene, nail polish remover, or spray).**

Failure to comply may cause an explosion or fire.

- **Do not aim the flash directly at your eyes.**

Failure to comply may cause injury to your eyes.

- **Ensure that neither you nor anyone else looks directly into the treatment window of the device.**

Failure to comply may cause injury to your eyes.

- **Do not use the device in a dark room.**

Failure to comply may cause injury to your eyes.

- **Do not direct the treatment window towards a black object.**

Failure to comply may result in smoke, fire, device malfunction, or injury.



Prohibited

- **Do not apply the flash on clothing, head hair, or objects.**

Failure to comply may cause scorching, discoloration, or damage.



Prohibited

- **Refrain from applying flash treatments without first removing foundation and other makeup products that have been applied to the skin.**

Be sure to remove makeup before using the device. Failure to comply may cause a burn injury or skin problems.

- **Do not use the device when the fan is not working (when a humming sound cannot be heard).**

Failure to comply may cause smoke, fire, malfunction, or injury. Please contact the store where you have purchased.

- **Do not pull out the plug from the outlet or from the device during use.**

Failure to comply may cause an electric shock or fire due to a short circuit.

- **Do not use or store the device in the following locations.**

- A high-humidity location such as a bathroom
- A dusty location
- A location where there is exposure to direct sunlight
- A location where water or other liquids may splash on the device
- A slope, locations subject to vibrations or shocks, outdoor locations
- A location where there are strong radio waves or magnetism
- A location where there is exposure to oil smoke or steam
- A location that becomes hot, such as near a heating device (for example, stoves)
- A location with drastic temperature changes (near an air conditioner or heater)

Failure to comply may lead to fire, electric shock, accident, malfunction, or device failure.

WARNING (Use and storage)

- Do not overuse the device. For usage time and frequency, follow the instructions in the User's Manual.
Failure to comply may cause physical injury.



Prohibited

- Do not leave the device unattended while it is turned on.
Doing so may cause the device to catch fire or may lead to an accident.

- Do not use the device for any purpose other than those stated in the User's Manual.

Failure to comply may cause an injury or a failure of the device.



Non-intended
uses are
prohibited.

- Use by unattended children or persons with physical disabilities should be avoided. A minor should use the device only under the supervision of his or her guardian who has given their consent. Do not allow children to play with the device.

Failure to comply may cause an accident or injury.



User
restrictions

- Do not use the device on people who are unable to express their will or who have reduced skin sensitivity.
Failure to comply may cause an accident or injury.
- Keep the device away from children and do not allow them to use the device.
Failure to comply may cause an accident or injury.

- Do not use the device in a bathroom or near a bathtub, shower, washbasin, or other vessel filled with water.

Failure to comply may cause an electric shock, short circuit, fire, or failure of the device.



Exposure
to water
is prohibited

- Do not splash water on the device (for example, the main unit) or immerse it in water. Do not use the device in a location where it may be splashed with water or in a high-humidity location.
Failure to comply may lead to fire, electric shock, or malfunction.

- Do not use the device outdoors or in direct sunlight.
Doing so may cause a failure of the device.



Outdoor use
is prohibited

- Do not place heavy items on the device, drop the device during use, or otherwise subject the device to a strong blow or force.

Failure to comply may cause an injury, a failure of the device, or the development of cracks in the glass tube.



Pumping
is prohibited

WARNING (Use and storage)

- **Do not disassemble, modify, or attempt to repair this device in any way.**

Failure to comply may cause a fire, electric shock or injury. Contact the store where you have purchased for advice on repairing the device.



Disassembly
is prohibited

- **Stop using the device and do not use it again when an odor or abnormal generation of heat is detected.**

Failure to comply may cause the device to become damaged or catch fire. Please contact the store where you have purchased.



Use is prohibited
when an odor
or abnormal
generation of
heat is detected

- **Do not leave the device near fire, in the hot sun, in a bathroom or other high-temperature/high-humidity location, or near combustible gas or a flammable object. Do not throw the device into a fire or bring it close to a heater.**

Failure to comply may result in the generation of heat, fire, the explosion of the device, or failure of the device.



Combustion
warning

- **Be sure to check the plug configuration and the voltage of the power outlet before using.**

Failure to comply may result in the generation of heat, fire, or a malfunction of the device.

- **Do not put the device into a microwave or pressure vessel.**

Failure to comply may result in the generation of heat, fire, or an explosion.



Burst

- **If water or foreign matter gets into the device (main unit, etc.), immediately stop using it and request an inspection.**

Failure to comply may lead to fire or electric shock.



Compliance
is mandatory

- **If you feel an earthquake, see lightning, or hear thunder, turn off the power and unplug the device.**
- Failure to comply may lead to fire or electric shock.

- **If the device is dropped or damaged, immediately turn off the power, pull out the plug, and be sure to request repair/inspection.**

Failure to comply may lead to fire or electric shock. Please contact the store where you have purchased.

- **Regularly check and clean the treatment window to make sure it is clean.**

Failure to comply may cause a burn injury or a failure of the device.

ATTENTION (Power supply)

- If a power failure occurs, turn off the power immediately and disconnect the AC adapter from the outlet.

Failure to comply may lead to an accident.



ATTENTION (Use and storage)

- Do not illuminate the same area more than once a day.

·Excessive application of flashes will not produce a better outcome. Please observe the "Treatment cycle" (page 18).
·Slight overlapping of flash application on areas during a single treatment does not pose a problem.

- Do not illuminate the same area for 5 or more minutes per day.

If you perceive a skin abnormality, immediately stop using the device. Consult a physician if the abnormality does not subside.

*One body part: (examples) upper arm, forearm, thigh, lower leg, underarm, face

- Before use, check to make sure that the treatment window is not cracked, deformed, or broken. Do not use the device if it is damaged or malfunctioning.

Failure to comply may cause your skin to suffer an injury or other problems, or lead to existing conditions worsening.

- Do not use perfume or cosmetics on the treated area immediately after treatment is completed. (excluding any cream or lotion used for moisture retention.)

Failure to comply may cause skin problems.

- Do not wipe the device with solvents such as alcohol, thinner, benzene, nail polish remover, kerosene, chemicals, or polishing powder.

Failure to comply may cause damage to the surface of the device.

- Do not apply an antiperspirant immediately after using the device on the underarm area.

Avoid applying an antiperspirant until any redness completely disappears.

- Do not use the device on pets.

Failure to comply may cause an injury.



ATTENTION (Use and storage)

- **Do not press the device strongly or use it intensively on a single spot.**
Pressing the device more strongly than necessary or using it intensively on a single spot may damage the skin.
- **Do not bring the device into contact with another beauty device or electronic appliance while the device is powered on.**
Do not use with other electric products (electric blankets, bed sheets, "kotatsu" foot warmers, etc.) at the same time.
Failure to comply may lead to an accident due to malfunction.
- **Do not use the device with your clothes on.**



-
- **When using the device on your face, be sure to use it in Face mode.**
Failure to comply may cause skin damage.
 - **Make sure to remove makeup before using the device on your face.**
If you use the device over makeup, the foundation may cause dark staining of the skin.
 - **Use only Level 1 upon first use and until you are accustomed to use.**
Even if you do not feel any pain, make only small adjustments to the intensity level. Use the device at an intensity level that suits you.
Make only small adjustments to the intensity level. Failure to comply may cause an accident or injury.
 - **Be sure to comply with the following items.**
Failure to comply may cause an accident, injury, or a failure of the device.
 - After sunbathing, do NOT use the device in combination with peeling, skin packs, or other treatments.
 - Be sure to turn the device off after use.
 - **Wear sunglasses if you feel that the light is too bright during your treatment.**
Failure to comply may cause injury to your eyes or skin.



ATTENTION (Use and storage)

- During treatment, carefully check that your skin is not showing any abnormalities. Discontinue use immediately if you start to feel pain or fatigue or if you see any skin reaction during use. Stop using the device immediately if the treatment does not agree with your skin. If symptoms continue the following day and beyond, consult a physician.
- A minor should safely and carefully use the device only under the supervision of his or her guardian, and immediately discontinue use whenever an abnormality arises.
Due to factors related to hormonal balance, damage may be caused to the skin, or the treatment may not yield observable results if the user has not reached full adulthood.
- Moisturize your skin daily and otherwise prevent your skin from drying out.
If your skin is dry, it is possible that the device may not work properly.
- Be sure to remove excess hair with a handheld shaver or razor before carrying out flash treatment.
Applying flash treatment without first removing excess hair may cause bits of hair to adhere to the device and produce scorch marks.
* Avoid tweezing or removing hair with wax or other such means.
- Be sure to remove any bits of hair remaining on the skin.
Your skin may feel hot when the flash is applied to hair debris.
- After each use, refer to 'After Treatment' (page 25) and remove bits of hair, cosmetics, and other particles that have adhered to the device.
Failure to comply may cause scorch marks, a failure of the device, or skin problems.
- Carefully clean the parts of the device that come into contact with the skin every time you use it, and keep the device in clean at all times.
Failure to comply may cause skin problems.



Compliance
is mandatory

ATTENTION (Use and storage)

- **When cleaning the device (main unit, etc.), be sure to turn off the power, disconnect the AC adapter from the outlet, and follow the instructions in the User's Manual.**

Failure to comply may lead to an accident or malfunctions.

- **Avoid exposing the treated area to ultraviolet rays for 24 hours after flash treatment.**

If redness continues to be exhibited even after 24 hours, avoid UV light. Also, protect your skin with sunscreen while undergoing treatment.

Failure to comply may cause skin problems.

- **If your skin remains hot after use, cool your skin down sufficiently by applying a cold towel.**

Your skin may feel hot after carrying out treatment depending on your physical constitution and method of use. If you are concerned about the hot flashes on your skin, cool it down by applying a cold pack or cold towel to the area in question. Consult a physician if the abnormality does not subside.

- **Ensure that the temperature of the usage environment is properly maintained.**

The range of environmental temperature in which the device can be used is from 5 to 35° C. Verify the ambient temperature of the location where the device will be used. The device may not correspond effectively to the ambient temperature of a new location immediately after being moved there. Failure to comply may cause the device to fail or malfunction.

- **Refrain from taking a bath (showering is permitted), exercising, or drinking alcoholic beverages on the day of treatment.**

Failure to comply may cause skin problems.

- **When you have not used the device for a long period of time, check to see that it is functioning normally before re-starting use.**

Failure to comply may lead to an accident or malfunctions.



Compliance
is mandatory

When Using Outside of Japan

WARNING

- Use only at AC100-240V to avoid a potential fire or electrical shock hazard.

Failure to comply may cause a fire or electric shock.



- Be sure to check the plug configuration and the voltage of the power outlet before use.


Failure to comply may result in the generation of heat, fire, or a malfunction of the device.



- The plug configuration and voltage vary depending on the region and the type of accommodation. Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a travel agency or a consumer electronics retailer.

Names and Functions of Components

Main Unit

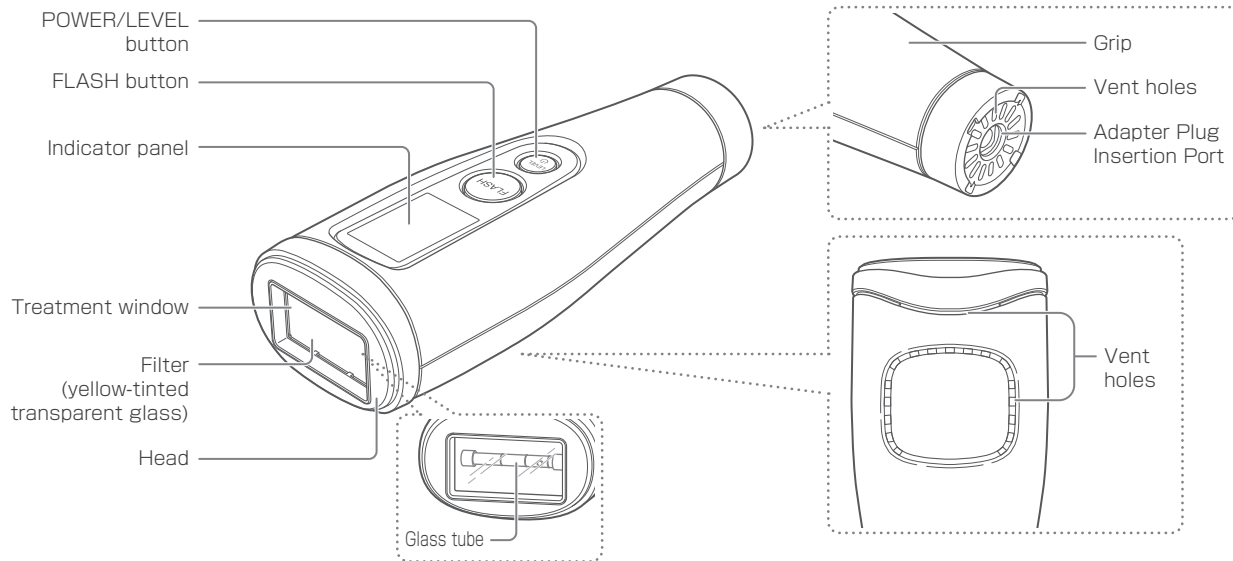


Accessories

● Face attachment

● AC adapter

● User's Manual

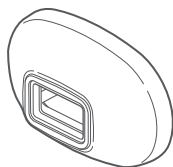

Use when performing treatment on the face (page 15).

The design and specifications of the device and peripheral accessories are subject to change for improvement purposes without notice.

Functions

Safety functions

● Touch sensor

You can flash your skin once the device detects that the treatment window is fully seated over your skin.

Auto Power OFF

About ten minutes after the device has been powered on, a buzzer sound will be emitted twice before the power is automatically shut off.

Cooling feature

After turning the power off, the fan will continue to operate for 90 seconds. This is a feature for cooling the device. Do not unplug the AC adapter until the fan stops.




· Always wait at least ten minutes after turning off the power before using the device again. Continual use may cause a failure of the device.



How to connect the face attachment

The face attachment can only connect to the device in one direction.

●Side view

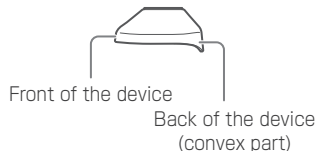

●Back (magnetic part)

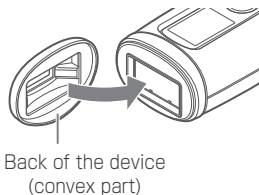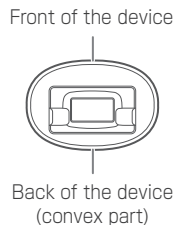


Check the direction of the face attachment and mount it from above so that it covers the treatment window on the device.

Since it is magnetic, it will stay in place if it is covering the treatment window.

* Make sure nothing is between the treatment window and the face attachment.

· If you mount or unmount the face attachment while the power is on, a buzzer will sound once and the mode will switch automatically.

How to remove the face attachment

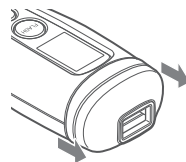


Hold the face attachment on both sides and pull it toward you.

Modes


You can perform skin treatment with the following modes using this device.

Body mode

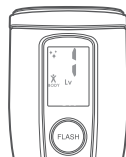

<Recommended areas of use>

●Arms

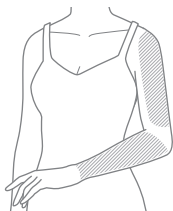

●Legs



●Underarms



●Bikini lines



Face Mode



The device will switch to Face mode automatically when the face attachment is mounted. Be sure to mount the face attachment and use Face mode when caring for the face.

<Recommended areas of use>

●Face (cheeks)



●Face (around the mouth)




●Fingers




* The device cannot be used on the lips.

Areas on which use of the device is prohibited

When performing treatment, do not use the device on the following areas.

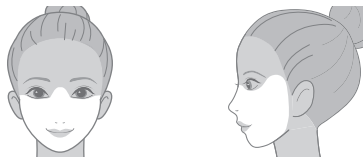

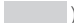- The head, scalp hair, nape of the neck, the forehead, brows, eyeballs, eyelids, around the eyes, ears, the Adam's apple, inside the mouth, lips (areas indicated in the above illustration as: )
- Nipples, areolae, inside the belly button, areas where a large vein or other blood vessels can be seen, sunburnt areas, the pubic region, or any mucous membranes.
- Acne, pigmentation or moles, injuries (scars or other dark-colored areas)
- Areas that have been subjected to a cosmetic procedure or surgery, tattoos or areas where the body has been painted, areas to which makeup has been applied.

Before Use

Before carrying out treatment, be sure to check the following contents.

Preparing for treatment

- Be sure to remove excess hair at the area to be treated with a handheld shaver or razor.

A burn injury may result if excess hair is allowed to grow long or remain on the surface of the skin at the area to be treated.

* Avoid tweezing or removing hair with wax or other such means.

- Be sure to perform a Flash test before using the device for the first time.

<Treatment test>



Apply a flash to an area near the intended care area starting with Level 1 in a single pulse. Treat the body using Body mode and the face using Face mode.

· You can use the device if no abnormality appears on your skin 24 hours after performing the flash test.

- Do not use the device while it is set to a level at which pain can be felt.
- Use the device while adjusting the level.

Since your perception will differ depending on the area where the device is used, adjust the level accordingly.

- Do not use the device while it is set to an unreasonably high level.

In particular, the skin of a darker-skinned person or a person who has gotten suntanned will more easily absorb light. For such persons, failure to comply may cause skin problems.


● **Check to make sure there is nothing unusual about the head, treatment window, or attachment. Do not use if damaged or deformed.**

- Be sure to clean the device if it gets dirty. Failure to comply may damage the device or cause burns, injuries, or skin problems.
- See "How To Clean" (page 25) for usage instructions.

● **Lighting during use**

- Using the device in a very dark environment may cause injury to your eyes.
- The recommended environment is to use [on a table in a room with ceiling lighting].
- Avoid use in the following conditions.
 - An indirectly lit room
 - A very dark room, such as one lit only by miniature bulbs
 - Dark places such as under a table

Treatment cycle

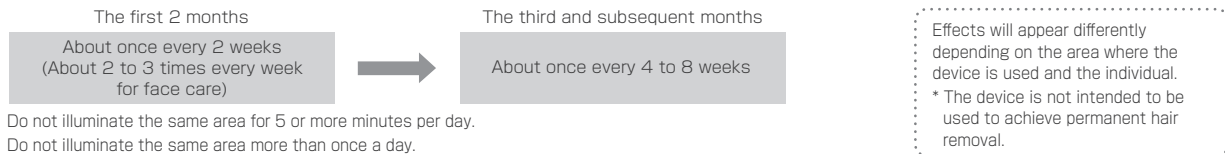

Do not illuminate the same area for 5 or more minutes per day.
Do not illuminate the same area more than once a day.



ATTENTION

● **Do not engage in excessive use of the device.**


If you wish to use the device in excess of the aforementioned frequency, be sure to use it while checking the condition of the skin.

In addition, always adhere to the prescribed frequency of use per day. Failure to comply may cause skin problems.

How to Use

Basic method of use

- Turn on the device**
Insert the AC adapter plug into the device and the AC adapter into a wall outlet and hold down the POWER/LEVEL button (approximately 3 seconds).
Standby mode is complete when the indicator panel is illuminated.

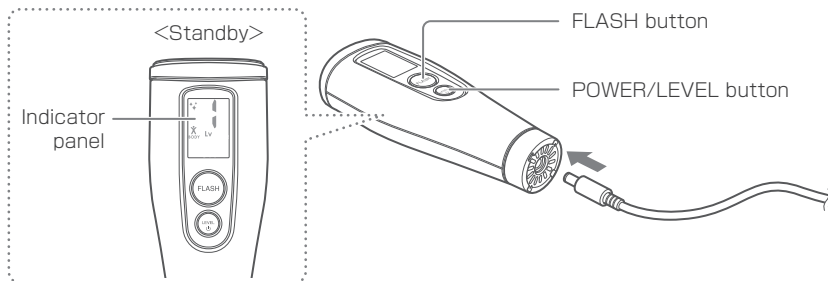

- The device will switch to Face mode when the face attachment is mounted.



- Do not overload an electrical circuit with the AC adapter.
Be sure that the AC adapter is inserted directly into a wall outlet.

2


Select the level

Every time you press POWER/LEVEL button, until the buzzer sounds once to change the level.

Body mode

Level switching will occur according to the following pattern.


→ Level 1 → 2 → 3 → 4 → 5
(When the power is turned on)

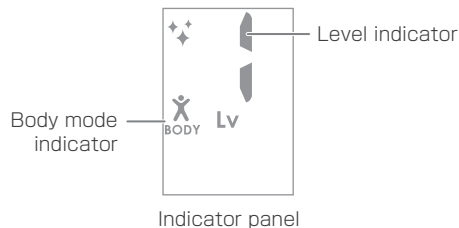

Face mode

Level switching will occur according to the following pattern.

→ Level 1 → 2 → 3 → 4 → 5
(When the power is turned on)

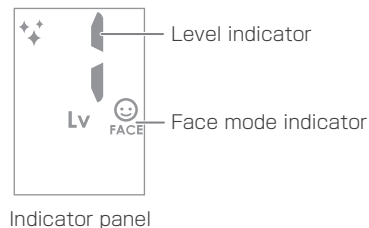

· Use only Level 1 upon first use and until you are accustomed to use.

3


Treatment

Press the entire treatment window against the skin. Preparation for treatment is complete when 2 short beeps sound and the OK lamp on the indicator panel is illuminated.

Press the FLASH button to use a flash.

See 'Points For Your Treatment' (page 24) for tips on performing treatments for different areas.

Both Body and Face modes

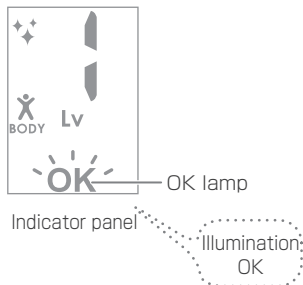

* See the same illustrations used for Body mode.

● Single pulse (one time)



Press the FLASH button to administer single pulses of illumination while shifting the position of the device.



- Each single pulse of illumination requires a few seconds of charge time.

● Continual treatment (up to 10 pulses)




* In Body mode only. Continual treatment is not possible in Face mode.

With the power on, hold the FLASH button down for about 3 seconds until it beeps twice to turn continual treatment on/off. Move the device over the skin with each pulse of illumination.

- * Turn the device on and off while the treatment window is removed from the skin.
- Continual treatment will end when the device beeps 3 times when removed from the skin for about 5 seconds during continual treatment.
- Each time continual treatment is used (up to 10 pulses), a long beep will sound, and the device will enter Standby mode. Press the FLASH button for a short time with the treatment window pressed against the skin to resume continual treatment (up to 10 pulses).
Standby: The mark  and the OK lamp will blink.
- During continual use, the mark  will display on the indicator panel.
- The level cannot be changed during continual treatment.



- Change the position of the device little by little during treatment to treat the entire area.
- Also, move the treatment area during continual use.
- After treating one area, remove the device from the skin and move it to the next area.
- The size of the treatment window is the area of skin that can be treated at one time.
- Do not treat the same area more than once a day.
- Since the skin around the face contour is susceptible to getting caught in the treatment window, be careful not to press too strongly while flashing.

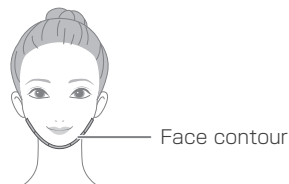

4

Turn off the device


- Hold down the POWER/LEVEL button (approximately 3 seconds).
- Verify that the fan has stopped (approximately 90 seconds) and disconnect the AC adapter from the device and the wall outlet.

After treatment, cool the skin using ice or a cool towel. Then, moisturize the area using a skin moisturizer.

Points For Your Treatment

Keep the following points in mind when performing treatment.

● Arms/legs

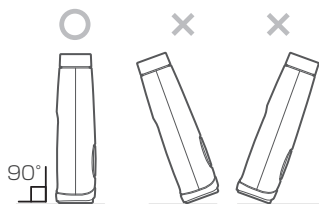

Hold the treatment window perpendicular to the skin.

● Underarms

Perform treatment while stretching out the skin with your fingers.

- Before performing treatment, carefully wipe away any antiperspirant or other cosmetic products on the skin.

● Bikini lines

Use the device at a low level on areas where there is pigmentation or areas with darker skin tone.

- * Be sure to start from Level 1 and adjust as you become more accustomed to the device.
- * The device cannot be used on external genitalia or surrounding mucous membranes.

● Face

Be sure to remove makeup before using the device.

- Use the device with a degree of intensity sufficient for you to perceive a slight warmth (ideally 2 to 3 times a week).
- When using the device on a narrow area such as mustache or beard, or a area with a sharp angle, perform the treatment while stretching out the skin with your fingers.
- Refrain from suntanning while undergoing treatment and avoid direct sunlight after treatment.

⚠ ATTENTION

● When performing treatment on the face, make sure to use Face mode.

Using Body mode could cause a burn injury.

● Ensure that neither you nor anyone else looks directly into the treatment window of the device.

Failure to comply may cause injury to your eyes.

● When treating the face, press the treatment window firmly against the skin.

If the treatment window is not pressed firmly against the skin, you may experience bright flashing.

After Treatment

How To Clean

Be sure to clean the treatment window after each use.

Regularly clean the main body as well to prevent the device from failing or malfunctioning.

WARNING

- **When cleaning the device, disconnect the AC adapter from the device and the wall outlet.**

Failure to comply may cause an electric shock, short circuit, or fire.

- **Clean the device in a manner that prevents water from entering it.**

Failure to comply may cause an electric shock, short circuit, or fire.

Treatment window

1

Use a soft piece of cloth or cotton swab to remove stains and bits of loose hair.

* Wipe away fingerprints and other stains.

Failure to comply may cause the filter to become scorched, a diminishment of flash performance, a failure of the device, or skin problems.

Main unit (including head and attachment)

1

Wet a piece of cloth with water and firmly wring it before wiping away stains from the surface of the device.

· If the device is very dirty, wipe its surface with a sponge or piece of cloth containing a diluted detergent solution.

Do not wipe the device with solvents such as alcohol, thinner, benzene, nail polish remover, kerosene, chemicals, or polishing powder.



Storing

- If you will not be using the device for an extended period, clean it before storing it (refer to 'How To Clean' on page 25).
- When the device is not in use, turn off the main body and disconnect the AC adapter from both the main body and wall outlet.
- Avoid using or storing the device near a stove or heater, in a bathroom or other high-temperature/high-humidity location, locations where there is exposure to direct sunlight, or unstable locations.
- Keep the device out of the reach of children.

Maintenance and Inspections

- If you have not used the device for an extended period, inspect it carefully and check that it works properly and safely before using it again.
- If the device is malfunctioning, go to 'Troubleshooting' (page 27) and attempt to identify the problem before arranging to have the device repaired. If the problem persists, contact the store where you have purchased with a proof of purchase.
- If the device malfunctions, appears damaged, or needs service, please contact the store where you have purchased.

Additional Information

Troubleshooting

If you perceive any abnormality while using the device, verify the following points before arranging to have the device repaired.


Problem	Possible Causes and Solutions
Does not turn on	→ Connect the AC adapter to the device and a wall outlet.
	→ Hold down the POWER/LEVEL button (approximately 3 seconds) to turn on the device.
The device cannot be operated.	→ Refer to 'Error Indications' (page 29).
The device shuts off	→ The device is set to shut off automatically (page 13). To resume using the device, wait at least 10 minutes before proceeding.
	→ Do not overload an electrical circuit with the AC adapter. Make sure that the AC adapter is fully inserted into a wall outlet.
No flash pulse is emitted	→ Preparations for flash emission have not been completed. Wait until 2 short beeps sound.
	→ Ensure that the treatment window is fully seated over the skin and wait for the OK lamp to turn on before removing your hair (page 21).
No flash pulse is emitted (During continual treatment)	→ Press the treatment window firmly on your skin (page 21).
	→ The device is in standby. Set the desired level with the POWER/LEVEL button (page 20).
No flash pulse is emitted. (During single pulse)	→ Press the FLASH button (page 21).
	→ The device is in standby. Set the desired level with the POWER/LEVEL button (page 20).

Problem	Possible Causes and Solutions
The flash pulse feels hot or painful.	→ Press the POWER/LEVEL button to lower the level (page 20).
	→ There is excessive hair. Make sure to get rid of excess hair with a handheld shaver or razor.
	→ Review the areas on which use of the device is prohibited (page 16).
	→ Your skin has entered the treatment window. Do not press your skin strongly.
Flash pulses are bright	→ The treatment window is not in full contact with the skin. Have the treatment window come into firm contact with the skin.
	→ Please wear sunglasses.
There is a buzzing sound.	→ This is the sound of the fan operating. This is not a malfunction or failure of the device.

Error Indications

If an error arises, it will be indicated in the indication panel on the device.

Check the contents of the error in question as outlined below and eliminate the cause of the error accordingly.

Error indication	Solutions
	3 short buzzer sounds are emitted. The head is hot, so turn off the power and allow the device to rest for about 30 minutes.
E and 1, 2, 3, or 4 light up alternately	20 short buzzer sounds are emitted. Please contact the store where you have purchased (page 32).
E and 5 light up alternately	Ten short buzzer sounds are emitted. Please contact the store where you have purchased (page 32).

Frequently Asked Questions

Q The flash pulses feel hot and I'm experiencing redness or pain. Is this normal?

A

Start using the device with Level 1.

Once you have grown comfortable with this level, gradually raise the level.

Perceiving a slight measure of heat from flash pulses is normal and does not pose any problem.

This problem can be attributed to the presence of long bits of excess hair, or to the treatment of suntanned skin with flash pulses.

Use the device on skin that has not been suntanned and after removing excess hair with a handheld shaver or razor.

If you experience redness or pain, discontinue using the device immediately and cool down and moisturize your skin.

If there is no amelioration, discontinue using the device and see a physician.

* Be sure to use the face attachment when caring for the face. Failure to comply may cause a burn injury.

Q How often can I use the device per day?

A

Usage should be limited to 5 minutes per day for each body part.

* One body part: (examples) upper arm, forearm, thigh, lower leg, underarm, face

Do not illuminate the same area more than once a day.

Intense heat-generated stimulation may give rise to skin problems.

Q Immediately after a flash pulse is discharged, I can smell a burning odor. Is this normal?

A

The burning odor comes from the heating of excess hair caused by the flash pulse.

If a burning odor is emitted by the device itself, immediately discontinue use, turn the device off, and disconnect the AC adapter.

Product Specifications

Input rating	DC12V 3.0A
Power consumption	Approximately 3.3 W: Standby (Approximately 39 W: Instantaneous maximum)
Outside dimensions	Approximately 68 mm × 47 mm × 178 mm (W×D×H)
Weight	Approximately 270 g
Materials	Main unit: ABS resin, nylon resin, acrylic resin, glass, Face attachment: ABS resin, aluminum, stainless steel (treatment window), nylon resin
Country of origin	Japan
Accessories	Face attachment One AC adapter One User's Manual

After-sales Service

Disposing of the device

Follow your local government rules for proper disposal of this product.

Contacts for overseas customers

For countries and regions outside Japan, please contact the store where you have purchased your product for any inquiry about it. Alternatively visit YA-MAN Global Official Website (<https://www.ya-man.com/en/>).

Memo

Memo

