
waterpik®

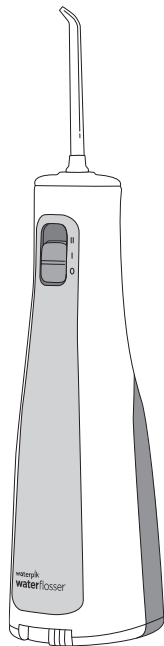**Waterpik® Water Flosser
WF-03**

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IMPORTANT SAFEGUARDS

Introduction

Thank you for your purchase from Ya-man LTD.

- Please read through this manual well to use this product safely.
- After reading this manual, please keep the manual to read it later whenever you need to do so.
- If you are giving this product to someone, please also give this manual too.

Caution for safety

- For safety please read this manual first before use.
- When you use the product, follow this manual.
- If you find something wrong with the product, please stop using it immediately and contact to our customer care for repair.
- Using defective product could be a risk to your oral health. Please do not use defective product.
- Information stated in this manual is very important for your safety. Please follow the instruction.
- To show the level of the hazard could be caused by faulty use of the product, we use these term "Danger" "Warning" and "Caution".

Hazard category

 **Danger:** imminent situation of people may decease or critically injured by faulty use.

 **Warning:** people may decease or critically injured by faulty use

 **Caution:** May cause injure to people or damage to their property.

*By property it means house, furniture, pets and domestic animals.

DANGER (use and storage of the product)

Prohibited

The following persons must not use this product.

There is the risk that the medical electronic device may malfunction. Be sure to consult with a physician before using this product.

- Persons who are using medical electronic devices
- Pacemakers and other implanted medical electronic devices
- Artificial heart-lung machines and other medical electronic devices used for life support
- Electrocardiographs and other electronic medical devices worn externally

WARNING (use and storage of the product)

Use only for the
designated
purpose

Do not use for purposes other than those specified in the instruction manual.

Always follow the instructions of your dental professional. There is the risk of accident or injury.

Must observe

Use only the tips and accessories recommended by Water Pik Inc.

There is the risk of accident or injury.

Must observe

The following people should consult with a dentist before using the product.

- If there is a suspicion of periodontal disease or you are being treated for symptoms of concern related to the teeth or oral cavity, contact a dentist.
- If you are taking antibiotics prescribed by a dentist, consult with the dentist before using this product or any other oral hygiene devices.

Prohibited

Do not use for other than the specified parts and the following parts.

Sensitive parts such as eyes, nose, ears, back of tongue, and throat. Since this product shoots water in strong pressure, there is a possibility of serious damage to these sensitive areas.

* If you cannot make your own judgment, consult with your doctor before use.

Must observe

Remove oral accessories before using this product.

There is the risk of injury.

Prohibited

Do not put hot water in the reservoir. Use room temperature water or warm water only.

There is the risk of burns.

Must observe

Use only the solutions recommended by a dentist.

There is the risk of injury and product failure.

Prohibited

**Do not use the product excessively and continuously.
(For business use etc.)**

There is the risk of heating and product failure.

Must observe

Observe the guideline of using the product twice a day.

Excessive use may cause damage to the tongue or oral cavity.

* This is only a guideline. Consider your personal health conditions carefully before using.

Prohibited

Do not use iodine, alkali salts (saline solution), or high-concentration essential oil that does not dissolve in water.

The performance of this product will be reduced, and the product life will be shortened. Also, it is not covered by the product warranty.

Fire caution

Do not place the product in fire, heat the product or use the product at high temperatures.

There is the risk of heating, ignition and rupture.

Prohibited

Do not immerse the product in water or leave it near water. Do not use the product or leave it in a bathroom or close to a vessel that contains water.

Water may enter inside and cause electric shock, short circuit, fire or malfunction.

Restriction
on use

Store the product in a location out of the reach of infants. Do not allow children to play with the product.

There is the risk of accident or injury.

IMPORTANT SAFEGUARDS

Restriction
on use

Do not use this product for people who cannot express their intentions or who have a weak oral sensation. Avoid operating the product only by people with physical difficulty. For minors, use under parental supervision and parental consent.

There is the risk of electric shock, injury and burns.

Prohibited

For reasons of sanitation, the product should not be used by more than one person.

There is the risk of injury and other problems.

Prohibited

Do not insert any foreign objects into the product opening.

There is the risk of product failure.

No disassembly

Do not disassemble, modify or repair the product.

There is the risk of fire, electric shock or injury. For repairs, contact the retailer where you purchased the product.

* This product does not contain any parts which can be repaired by the customer.

No use when
there is
abnormal odor

Do not use the product if it has unusual smell or overheating.

It may cause damage or fire.

No outdoor use

Do not use the product outdoors.

Do not operate the product in areas where aerosol (spray) products or oxygen inhaler are used.

There is the risk of injury and product failure.

Fire caution

Do not use or leave the product close to fire, exposed to direct sunlight or in a dusty location.

There is the risk of heating, ignition, rupture and malfunction.

No impacts

Do not place a heavy object on the product or drop it during use.

There is the risk of injury and product failure.

Do not drop

Be sure to place the product in a flat, stable location. Do not place it anywhere that is not stable.

There is the risk of injury, ignition and product failure.

Fire caution

Do not put the product in a microwave oven or high-pressure container.

There is the risk of heating, ignition and rupture.

Fire caution

When/After using the product, be sure to wipe off any moisture from the product using a dry cloth or similar means. Then, store it.

There is the risk of electric shock, short circuit, ignition and product failure.

Caution (use and storage of the product)

Prohibited

Do not use if this product is cracked, damaged or defective.

It may cause damage to the oral cavity.

Prohibited

Do not press the tip forcefully against your teeth or gums.

There is the risk of injury.

Prohibited

Do not wipe the product with solvents such as alcohol, thinner, benzene, or nail polish remover.

This may cause damage to the surface of the main unit.

Prohibited

Do not expose the product to hot water or hot air from a hair dryer.

It may cause the product to deform or melt.

Prohibited

Do not set the product down when the power is ON. Be sure to turn the power OFF after use.

There is the risk of injury and damage to furniture or other items.

Prohibited

Do not let the product contact with other electrical products while the power is on. Also, keep it away from heat sources such as electric heating equipment (electric blankets, kotatsu).

There is the risk of injury and product failure.

Prohibited

Do not use the product on pets.

There is the risk of injury.

Prohibited

Do not use the product if you have a cut or other abnormality inside your oral cavity.

There is the risk of injury.

Must observe

If you feel unwell during use or experience any abnormalities in your body, stop using it immediately.

If the problem continues after the next day, consult a doctor.

Must observe

After each use, please clean the product. Please refer to "Caring for the product" (P.12).

Without proper care for the product, it may cause oral troubles and breakdowns.

Cautions (Concerning the Battery)

Must observe

Be sure to observe the following precautions.

- This product uses AA alkaline dry-cell batteries or NiMH (nickel metal hydride) rechargeable batteries.
- Be sure to observe the precautions printed on the batteries.
- Check the orientation of the + and - poles before inserting the batteries.
- When replacing a battery, be sure to use a new battery of the same type and the same brand.
- Remove dead batteries promptly.
- Use batteries before their recommended expiration dates.
- When using NiMH (nickel metal hydride) rechargeable batteries, read the instruction manuals of the batteries and charger, and use them correctly.



There is the risk of injury or damage to surrounding property caused by battery heating, rupture, or fluid leakage.

The use of lithium batteries will void the product warranty. Such batteries may also damage the product.

Product Description	7
Getting Started	8
Cleaning Your Waterpik® Water Flosser	12
Troubleshooting	14
FAQ	15
Warranty Information	16

[Main unit]

1. Power/pressure switch
2. Power handle
3. Battery cover with lock mechanism
4. Tip
5. Tip eject button
6. Water filling port (cap)
7. Reservoir
8. Water intake port

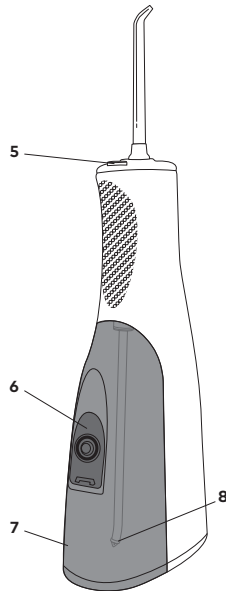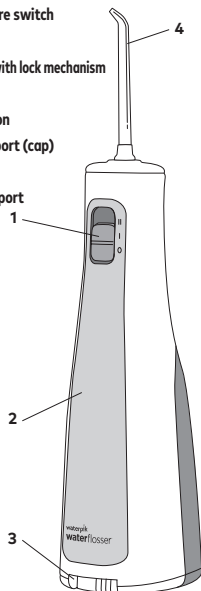4. The tips are accessories and are not covered by the warranty.

Tip types and oral care uses

A. Orthodontic Tip*

- Braces/General Use

B. Pik Pocket™ Tip*

- Periodontal Pockets/Furcations

C. Plaque Seeker™ Tip*

- Implants/Crowns/Bridges/Retainers/General Use

D. Toothbrush Tip*

- General Use



E. Classic Jet Tip**

- General Use

F. Tongue Cleaner**

- Fresher Breath

Not all of these tips are provided with all models.

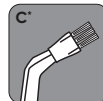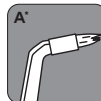*Replace Every 3 months **Replace Every 6 months

For details about the tips and other accessories, and about how to use the product, please see our company's website (<https://www.ya-man.com>).

Inserting Batteries

1. Slide the reservoir as shown in the figure to remove it.
2. Rotate the lock mechanism on the battery cover counterclockwise to the unlock position, then remove the battery cover from the main unit.
3. Insert AA batteries as shown in the figure. AA alkaline dry-cell batteries are provided with the product for trial purposes.
4. Press the battery cover in place so the arrow on the cover matches the arrow on the main unit, then rotate the battery cover lock knob clockwise.







If the water pressure drops while the product is operating, replace the batteries. The product can be used around 40 times after the batteries are replaced. However this varies depending on the battery type and conditions of use.

Important

When replacing the batteries, apply commercially available lubricant to the O-ring and battery case.

Apply a thin, even coat of commercially available lubricant to the O-ring (black waterproof rubber ring) on the inside of the battery cover.

Apply a thin, even coat of commercially available lubricant to the parts of the battery case cover on the main unit side which the O-ring contacts as shown in the figure.

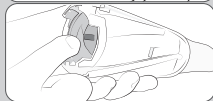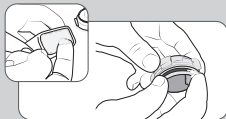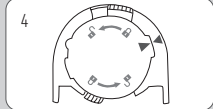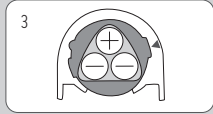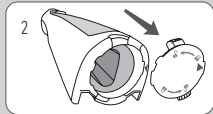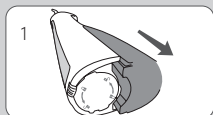It is recommended that this product be used after you brush your teeth. Guideline for use: 2 times/day

Filling the Reservoir

Open the reservoir cap and fill the reservoir with cold water or warm water.

Inserting and Removing Tips

Insert the tip into the insertion hole located on the top of the Water Flosser main unit. When the tip is locked correctly in the designated position, the colored part at the center of the tip will be perfectly aligned with the top of the main unit. To remove the tip from the handle, press the tip eject button and pull tip from handle.

Adjusting the Pressure Setting

Simply slide the power/water pressure adjustment switch to adjust the Water Flosser pressure setting.

Slide the power switch up and align it with (I) for low pressure. Align it with (II) for a high-pressure water jet. Use the product in accordance with instructions from a dentist.

Recommended Technique

Lean low over sink and place the tip in mouth. Aim the tip toward teeth, and turn the unit ON (I).




For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

Caution

Do not press the tip forcefully against your teeth or gums.

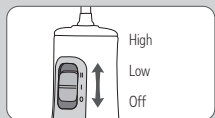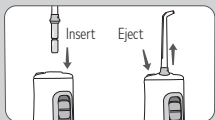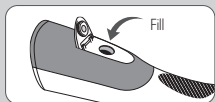
When Finished

When finished with use, slide the power/water pressure adjustment switch to "0" to turn the power OFF.

In order to keep the main unit and reservoir clean, empty the remaining liquid in the reservoir after each use. A small amount of water will remain inside the reservoir but this does not represent a problem for health.

Finally wipe off any water droplets on the main unit and reservoir using a towel or other dry cloth.

To install the reservoir back onto the main unit, slide the reservoir toward the top of the main unit.

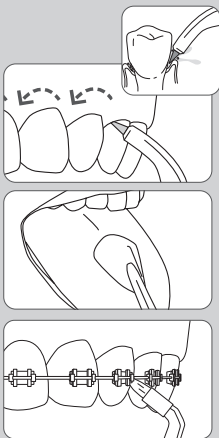TIP USAGE

Pik Pocket™ tip

The Pik Pocket™ tip is specially designed to direct a jet of water or solution deep into the pocket around the tooth.

When using the Pik Pocket™ tip, **set the main unit water pressure to the lowest setting** and angle the tip at 45 degrees from the tooth. Slowly insert the tip into the pocket around the tooth.

Tongue Cleaner tip

Set the water pressure to the lowest setting and place the tip near the center of your tongue. Gently press the tip against your tongue and pull it toward you. Adjust the water pressure to suit your preferences.


Orthodontic tip

When using the orthodontic tip, move the tip slowly and smoothly along your gums, then stop it for a moment to gently clean the space between the teeth and around the orthodontic bracket. Then move on to the next tooth.

**Toothbrush tip**

Insert the toothbrush tip into your mouth and place the brush head against your gums. It is not necessary to put toothpaste onto the toothbrush tip.

Use a gentle force (that does not bend the bristles) and move the brush back and forth and use short strokes like a manual toothbrush.

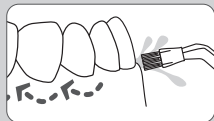**Plaque Seeker™ tip**

When using the Plaque Seeker™ tip, position the tip close to your teeth so that the bristles gently contact the teeth.

Move the tip slowly and smoothly along your gums, stopping it briefly at the spaces between the teeth. Brush gently and direct the water stream between the teeth.

**Caution**

Not all of these tips are provided with all models.

Using mouthwash or other solution

The Waterpik® Water Flosser can use mouthwash or other solution. In order to prevent clogging, after using a solution other than water, fill the reservoir halfway with warm water, then direct the tip toward the sink and operate the main unit until the reservoir is empty.

Caring for the product (daily)

● Main unit

Turn the power OFF, then use a soft cloth and a neutral detergent that does not contain a polishing agent, and wipe the product as necessary.

● Reservoir

Pour out the water in the reservoir and dry it in a well ventilated location. If there is dirt on the surface, wipe with a soft cloth and a neutral detergent that does not contain a polishing agent.

The reservoir can be washed in a dishwasher if it is placed on the top rack. (The tips cannot be washed in a dishwasher.)

● Tips

Wipe with a soft cloth and a neutral detergent that does not contain a polishing agent.

Caution

- Do not wipe the main unit with alcohol, paint thinner, benzene, nail polish remover or other solvent. These may damage the product surface.
- Do not immerse the main unit in water.

Cleaning inside the main unit (once every 1 - 3 months)

Depending on the mineral content of the water that is used, water scale may build up inside the main unit.

If it is not removed, it may decrease product performance.

1. Fill the reservoir with water or warm water and add 2 tablespoons of white vinegar or grain vinegar.
2. Direct the handle and tip toward the sink and operate the main unit until the reservoir is empty.
3. Fill the reservoir with clean water or warm water and repeat the same procedure around 3 times to rinse out the main unit.

Storing the product

- After use, be sure to fully wipe off any water and dry the product before storing it. If the product is left with water on it, it may cause corrosion inside the main unit.
- If the product will not be used for a long period of time, clean the main unit before storing it.
- If the product will not be used for a prolonged period, remove the AA batteries.
- Do not place the product close to fire, exposed to direct sunlight, in a bathroom or other location exposed to high temperatures and humidity or in an unstable location. Store it in a clean location with good ventilation.
- Store the product in a location out of the reach of children.

Maintenance and inspection

- If the product was not used for a long period of time, be sure to check that the product operates correctly and safely before using it.
- If there is any problem with operation of the product, before requesting repairs, refer to the "Troubleshooting" (P.14) section and check the trouble. If the problem is still not corrected, contact the retailer where you purchased the product, or prepare the purchase invoice and contact the YA-MAN Call Center (P.16).
- If the product has malfunctioned or is damaged, contact the retailer where you purchased the product.

Handling and Disposal of AA batteries

- This product uses alkaline dry-cell batteries or NiMH (nickel metal hydride) rechargeable batteries.
- Check the orientation of the + and - poles and insert the batteries facing in the correct direction.
- Remove dead batteries from the Water Flosser and dispose of them in accordance with local waste sorting rules.
- The supply terminals are not to be short-circuited.
- Do not dispose of batteries in fire.
- The type and quality of the AA batteries used will affect product performance.
- When replacing a battery, be sure to use a new battery of the same type and the same brand.
- Replace batteries when flosser is no longer effective.
- Remove batteries when the appliance will not be used for a long period of time.
- The lifetime of ordinary batteries is around 40 uses.
- Be sure to observe the precautions printed on the batteries.
- Remove dead batteries promptly.
- Use dry-cell batteries before their recommended expiration dates.
- When using NiMH (nickel metal hydride) rechargeable batteries, read the instruction manuals of the batteries and charger, and use them correctly.

Caution

The use of lithium batteries will void the product warranty. Such batteries may also damage the product.

Symptom	Possible Causes	Solutions
Unit won't turn on?	<p>Typically occurs when unit goes into "lock out mode". Lock out mode can happen when:</p> <ul style="list-style-type: none"> - AA batteries have not been inserted or are dead. - the tip on the unit is clogged - the temp gets lower than 30-40 degrees F - The water intake port (on the end of the white water intake hose) inside the main unit is clogged. The device was dropped or was left unused for long period of time. 	<ul style="list-style-type: none"> - Replace with new batteries. - Operate with the tip removed. If it operates, replace the tip. - Operate with the reservoir removed. If it operates, perform cleaning of the water intake port (on the end of the white water intake hose).
Low or No water pressure	The tip is clogged, water intake port (on the end of the white water intake hose) is clogged, an additive was used, cleaning was insufficient or preparation for use was insufficient.	<ul style="list-style-type: none"> - Perform cleaning.
Tip Broken in Unit	The product was dropped.	<ul style="list-style-type: none"> - Screw a metal screw into the end of the broken tip. - Hit eject button and remove screw. Tip will come out with the screw.
Unit is leaking	The tip is clogged, the reservoir drain port is clogged, an additive was used, cleaning was insufficient or preparation for use was insufficient.	<ul style="list-style-type: none"> - If there is visible cracking of the reservoir, replace the reservoir. Contact the YA-MAN Call Center (P.16).
Tip Ejects while in use	Most commonly caused by build-up, defective tip, or defective unit.	<ul style="list-style-type: none"> - Try a different tip to see if problem still occurs - If you used additives like mouthwash or peroxide, you will want to flush water in the hole where tip is inserted.

FAQ

Where can I buy replacement tips and accessories?

Contact the YA-MAN Call Center. (P.16).

What is the warranty on this product?

1 Year

How long do the batteries last?

When new alkaline dry-cell batteries are inserted, they can be used around 40 times.

Can I use solutions in my water flosser?

Mouthwash or other solution can be used. Each time after such a solution was used, fill the reservoir with clean water and perform cleaning of the device. Failure to rinse the unit may result in loss of product performance and/or function. Growth of bacteria could also be enhanced. Note: iodine, saline or water insoluble concentrated essential oils cannot be used.

How long does the water last? How much does the reservoir hold?

The reservoir capacity is approximately 148 mL (enough for approximately 30 seconds of use).

What is the pressure range?

45-75 PSI with 2 pressure settings

WARRANTY INFORMATION

Disposing

- When disposing of this product, please dispose of it according to the rules of your area.

Product specifications

Power supply type	Dry-cell battery type
Dry-cell batteries	AA dry-cell batteries x 3 * Alkali or nickel metal hydrogen
External dimensions	Approximately W64 x D92 x H299 (mm) when tip is installed
Product weight	Approximately 340g
Reservoir capacity	Approximately 148 mL (enough for approximately 30 seconds of use)
Accessories	Jet tips: 2, orthodontic tip: 1, AA alkaline dry-cell batteries (for trial): 3, instruction manual (with warranty certificate)
Material	Main unit: ABS plastic, reservoir: copolyester, tips: copolyester/ABS plastic
Country of manufacture	China

Contact Information

If you have any questions about the product warranty, repairs or other matter, please contact the Call Center below. Or visit www.ya-man.com for more information.

YA-MAN Call Center

Hours: 9:00 - 17:00 (closed on Saturdays, Sundays, holidays, and during the year-end and New Year holidays)

From a land line or similar phone

 **0120-776-282** (Call is free.)

From a mobile phone or PHS

 **0570-550-637** (A call fee is charged.)

* When requesting repairs, prepare the warranty certificate or purchase invoice that contains the date of purchase and retailer name.

* The Call Center lines are often busy on Mondays.

Warranty information

A warranty for 1 year from the date of purchase is provided and covers problems that occur during ordinary use. (It does not cover consumable supplies or accessories.)

- After the warranty period expires, product shipping expenses will be borne by the customer.
- In the case of problems occurring after the warranty period expires, the customer will be charged for the cost (actual expense) of repairs.
- There are cases when a charge for repairs will be incurred even during the warranty period. Please read the explanation of the warranty terms carefully.
- In the event of a problem with the product, contact the retailer or the YA-MAN Call Center (P.16).

Warranty Terms

1. Problems which occur under conditions of normal use as prescribed in the instruction manual, labels attached to the product, and other precautions shall be repaired at no charge by the retailer where the product was purchased.
2. To receive free repairs of a problem that occurred during the warranty period, present the product and this warranty certificate, and request repairs by the retailer where the product was purchased.
3. If your address will be changed, please consult in advance with the retailer where the product was purchased.
4. If the product was a gift and it is not possible to request repairs by the retailer where the product was purchased listed in this certificate, please contact the YA-MAN Call Center.
5. Even during the warranty period, a charge shall be required for repairs in the following circumstances.
 - (a) The problem or damage is the result of incorrect use or unauthorized repair or modification.
 - (b) The problem or damage is the result of the product tipping over or falling after it was purchased.
 - (c) The problem or damage is the result of fire, earthquake, flood or other disaster, pollution or abnormal voltage.
 - (d) The problem or damage is the result of an abnormality in another connected device.
 - (e) Parts such as batteries, consumable supplies or accessories require replacement.
 - (f) The problem or damage occurred when the product was used in a location not an ordinary household (for example, prolonged use for commercial purposes).
 - (g) The problem or damage occurred when the product was mounted on a vehicle, ship or similar location.
 - (h) The problem or damage is the result of product transport, relocation, falling or similar cause after it was purchased.
 - (i) The warranty certificate was not presented.
 - (j) The date of purchase or the retailer name or stamp is not present in the warranty certificate, or the contents of the warranty certificate were modified.
6. This warranty is valid only within Japan.

* This certificate is a pledge to provide free repairs based on the warranty period and terms printed on it. Therefore this warranty certificate does not limit the rights of the customer under the law. If you have any questions regarding repairs after the end of the warranty period or other matters, please contact the retailer where the product was purchased.

* A charge will be incurred for repairs in cases of deterioration or consumption resulting from use, or for replacement of consumable supplies.

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©2020 Water Pik, Inc.
Water Pik, Inc. a subsidiary