
waterpik®

Waterpik® Water Flosser WP-600 Series

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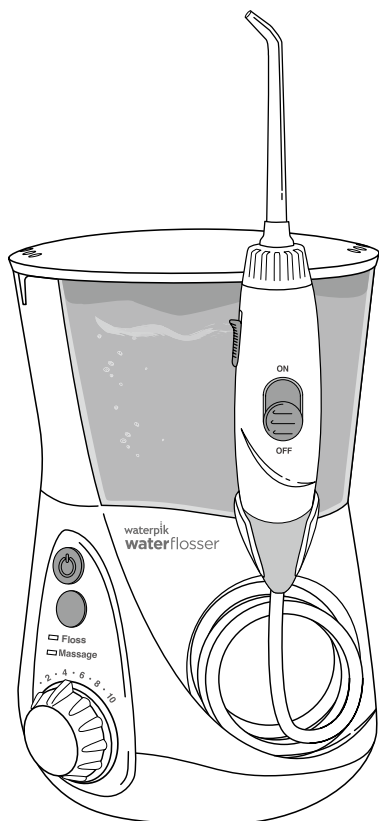IMPORTANT SAFEGUARDS

INTRODUCTION

Thank you for your purchase from Ya-man LTD.

- Please read through this manual well to use this product safely.
- After reading this manual, please keep the manual to read it later whenever you need to do so.
- If you are giving this product to someone, please also give this manual too.

CAUTION FOR SAFETY

- For safety please read this manual first before use.
- When you use the product, follow this manual.
- If you find something wrong with the product, please stop using it immediately and contact to our customer are for repair.
- Using defective product could be a risk to your oral health. Please do not use defective product.
- Information stated in this manual is very important for your safety. Please follow the instruction.
- To show the level of the hazard could be caused by faulty use of the product, we use these term "Danger" "Warning" and "Caution".

HAZARD CATEGORY

 **Danger:** imminent situation of people may decease or critically injured by faulty use.

 **Warning:** people may decease or critically injured by faulty use

 **Caution:** May cause injure to people or damage to their property.

*By property it means house, furniture, pets and domestic animals.

DANGER

Prohibited

The following persons must not use this product.

There is the risk that the medical electronic device may malfunction. Be sure to consult with a physician before using this product.

- ◆ Persons who are using medical electronic devices
 - Pacemakers and other implanted medical electronic devices
 - Artificial heart-lung machines and other medical electronic devices used for life support
 - Electrocardiographs and other electronic medical devices worn externally

WARNING (Power supply)

Must observe

Use the product in 100-240v. Free repeating warranty is only valid when products are used in Japan.

If not, it may cause fire, electric shock and burn.

Must observe

Plug the AC adaptor fully before use.

If not, it may cause fire, electric shock and fire.

Unplug

Hold power plug when you disconnect it. Do not grab the cable to pull out. Please make sure your hands are dry when you do so.

It may cause fire, electric shock, burn or cable disconnection.

Unplug

After using the product, be sure to turn the power OFF and unplug the power plug from the outlet.

There is the risk of fire, ignition and burns.

Must observe

Unplug the power plug when you clean the product.

There is the risk of electric shock, fire, ignition and burns.

Prohibited

Do not drag the power cable to move or lift the product.

It may cause electric shock, ignition and cable disconnection

Prohibited

Do not use the product with its cable bundled.

The heat from the cable will not be released and may cause fire.

Prohibited

Do not touch the power plug with wet hands.

There is the risk of electric shock, injury and burns.

Do not use a damaged cord or plug

Do not use when Power plug or cable is damaged.

There is the risk of electric shock, short circuit, ignition and product failure. If the power plug or cord is damaged, be sure to request repairs from the retailer where you purchased the product.

Unplug

In the event of any abnormality, turn OFF the product power and immediately unplug the power plug.

Do not use the product when there are any abnormal conditions such as smoke, abnormal noise or abnormal odor.

Prohibited

Do not allow the power plug and power cable to become wet (do not wash them or immerse them in water). Make sure the power plug is dry before plugging in the product.

There is the risk of electric shock, short circuit, ignition and product failure.

Prohibited

Do not damage, break, modify, overly bend, pull, twist or bundle the power cable. Do not put a heavy object on the cable. Do not allow the cable to become pinched in other objects.

There is the risk of breakage.

Must observe

When earthquake or lightning happens, turn off the product and unplug the unit.

Piggybacked outlets are prohibited

Do not overload an electrical circuit. It may cause electric shock and fire.

There is the risk of electric shock, short circuit and ignition.

WARNING (Use and storage)

Use only for the designated purpose

Do not use for purposes other than those specified in the instruction manual.

Always follow the instructions of your dental professional. If not, it may cause an accident or injury.

Must observe

Use only the tips and accessories recommended by Water Pik Inc.

Other items may cause accident or injury.

Must observe

The following people should consult with a dentist before using the product.

- If there is a suspicion of periodontal disease or you are being treated for symptoms of concern related to the teeth or oral cavity, contact a dentist.
- If you are taking antibiotics prescribed by a dentist, consult with the dentist before using this product or any other oral hygiene devices.

Prohibited

Do not use for other than the specified parts and the following parts.

Sensitive parts such as eyes, nose, ears, back of tongue, and throat. Since this product shoots water in strong pressure, there is a possibility of serious damage to these sensitive areas.

* If you cannot make your own judgment, consult with your doctor before use.

Must observe

Remove oral accessories before using this product.

It may cause injury.

Prohibited

Do not put hot water in the reservoir. Use room temperature water or warm water only.

There is the risk of burns.

Must observe

Use only the solutions recommended by a dentist.

There is the risk of injury and product failure.

Prohibited

Do not use the product excessively and continuously. (For business use etc.)

There is the risk of heating and product failure.

Must observe

Observe the guideline of using the product twice a day.

Excessive use may cause damage to the tongue or oral cavity.

* This is only a guideline. Consider your personal health conditions carefully before using.

Prohibited

Do not use iodine, alkali salts (saline solution), or high-concentration essential oil that does not dissolve in water.

The performance of this product will be reduced, and the product life will be shortened. Also, it is not covered by the product warranty.

Fire caution

Do not place the product in fire, heat the product or use the product at high temperatures.

There is the risk of heating, ignition and rupture.

Prohibited

Do not immerse the product in water or leave it near water. Do not use the product or leave it in a bathroom or close to a vessel that contains water.

Water may enter inside and cause electric shock, short circuit, fire or malfunction.

Restriction on use

Store the product in a location out of the reach of infants. Do not allow children to play with the product.

There is the risk of accident or injury.

Restriction on use

Do not use this product for people who cannot express their intentions or who have a weak oral sensation. Avoid operating the product only by people with physical difficulty. For minors, use under parental supervision and parental consent.

If not followed, it may cause electric shock, injury or burns.

Prohibited

For reasons of sanitation, the product should not be used by more than one person.

There is the risk of injury and other problems.

Prohibited

Do not insert any foreign objects into the product opening.

There is the risk of product failure.

No disassembly

Do not disassemble, modify or repair the product.

There is the risk of fire, electric shock or injury. For repairs, contact the retailer where you purchased the product.

* This product does not contain any parts which can be repaired by the customer.

No use when there is abnormal odor

Do not use the product if it has unusual smell or overheating.

It may cause damage or fire.

No outdoor use

Do not use the product outdoors. Do not operate the product in areas where aerosol (spray) products or oxygen inhaler are used.

There is the risk of injury and product failure.

Fire caution

Do not use or leave the product close to fire, exposed to direct sunlight or in a dusty location.

There is the risk of heating, ignition, rupture and malfunction.

No impacts

Do not place a heavy object on the product or drop it during use.

It may cause injury or malfunction.

Do not drop

Be sure to place the product in a flat, stable location. Do not place it anywhere that is not stable.

There is the risk of injury, ignition and product failure.

Fire caution

Do not put the product in a microwave oven or high-pressure container.

It may cause overheating, ignition or explosion.

Fire caution

When/After using the product, be sure to wipe off any moisture from the product using a dry cloth or similar means. Then, store it.

If not, it may cause electric shock, short circuit, fire or malfunction.

CAUTION (Use and storage)

Prohibited

Do not use if this product is cracked, damaged or defective.

It may cause damage to the oral cavity.

Prohibited

Do not wipe the product with solvents such as alcohol, thinner, benzene, or nail polish remover.

This may cause damage to the surface of the main unit.

Prohibited

Do not expose the product to hot water or hot air from a hair dryer.

It may cause the product to deform or melt.

Prohibited

Do not set the product down when the power is ON. Be sure to turn the power OFF after use.

There is the risk of injury and damage to furniture or other items.

Prohibited

Do not let the product contact with other electrical products while the power is on. Also, keep it away from heat sources such as electric heating equipment (electric blankets, kotatsu).

There is the risk of injury and product failure.

Prohibited

Do not use the product on pets.

It may cause injury.

Prohibited

Do not use the product if you have a cut or other abnormality inside your oral cavity.

There is the risk of injury.

Must observe

If you feel unwell during use or experience any abnormalities in your body, stop using it immediately.

Pay special attention to your health condition when using this product. If the problem continues after the next day, consult a doctor.

Must observe

After each use, please clean the product. Please refer to "Caring for the product" (P.12).

Without proper care for the product, it may cause oral troubles and breakdowns.

Prohibited

Do not press the tip forcefully against your teeth or gums.


There is the risk of injury.

TABLE OF CONTENTS

Product Description	7
Getting Started	8
Cleaning Your Waterpik® Water Flosser	11
Trouble Shooting	14
FAQ	14
Limited Two-Year Warranty	15

[Main unit]

1. Dual purpose covered lid
2. Reservoir
3. Tip storage under lid
4. Power button
5. Mode selector button
6. LED mode indicators
7. Pressure control dial
8. Exchangeable tips
9. Rotating control knob
10. Tip eject button
11. Water ON/OFF switch
12. Handle
13. Water hose
14. Base

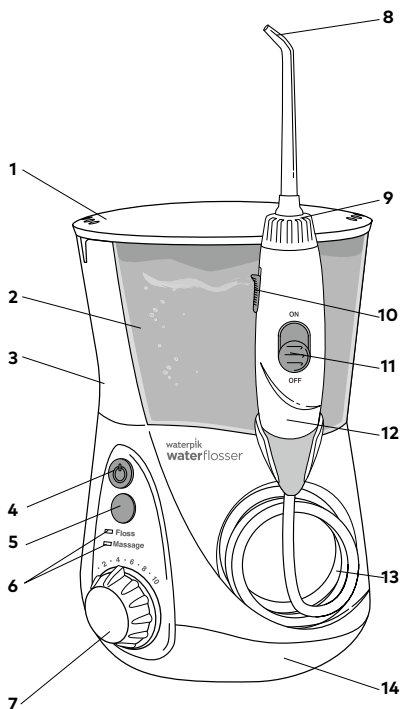8. The tips are accessories and are not covered by the warranty.

Tip types and oral care uses

A. Orthodontic Tip*

- Braces/General Use

B. Pik Pocket™ Tip*

- Periodontal Pockets/Furcations

C. Plaque Seeker™ Tip*

- Implants/Crowns/Bridges/
Retainers/General Use

D. Toothbrush Tip*

- General Use





E. Classic Jet Tip**

- General Use

F. Tongue Cleaner**

- Fresher Breath

NOT ALL TIPS INCLUDED IN ALL MODELS

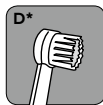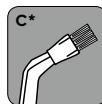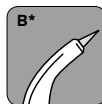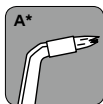*Replace Every 3 months **Replace Every 6 months

For details about the tips and other accessories, and about how to use the product, please see our company's website (<https://www.ya-man.com>).

Before using the product

It is recommended that this product be used after you brush your teeth.

Guideline for use: 2 times/day

Powering the Unit

Plug the cord into an electrical outlet.

Preparing the Reservoir

Lift up the reservoir lid, remove the reservoir, and fill it with room-temperature water or warm water. Return the reservoir to the main unit and press it down securely, then close the lid.

Inserting and Removing Tips





Insert tip into the center of the Water Flosser handle, pressing down firmly until the tip clicks into place. **When the tip is locked correctly in the designated position, the tip rotating knob will be perfectly aligned with the top of the handle.** To remove tip from the handle push the tip eject button on the handle.



Prime Unit (First Time Use Only)

Set pressure control dial on base to pressure setting No.10, point the handle and tip into sink. Turn ON the power button and turn the water stream ON/OFF switch to ON. Check that water flows out from the tip. If there is no problem, turn the switch OFF.

Adjusting the Pressure Setting

Turn the pressure control dial on the base of the unit to the lowest setting for first-time use. Gradually increase pressure over time to the setting you prefer. Use the product in accordance with instructions from a dentist.

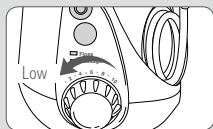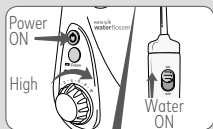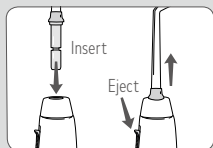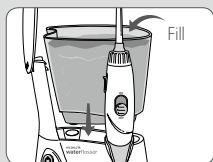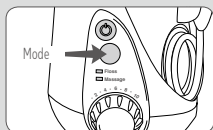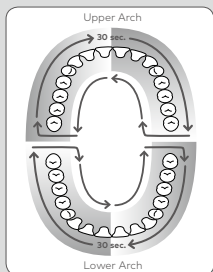Using the Water Flosser

Press the power button **and lean your body over the sink. Slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn water switch ON. When finished, turn water switch and power OFF.** During use you can stop the flow of water by sliding the water switch on handle to the OFF position.

⚠ Caution

- Be sure to place the tip in your mouth before turning the water stream ON/OFF switch to ON.
- Do not press the tip forcefully against your teeth or gums.

It is recommended that you start with your molars and work your way toward the front teeth. Insert the tip into your mouth and aim the water so that it contacts your gums perpendicularly. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

One-Minute Timer/30-second Pacer (Floss Mode Only)

Around 30 seconds after starting use, the product stops temporarily. This notifies you that you should move to the other arch of the oral cavity. After another 30 seconds (around 1 minute after the start of use), it stops temporarily again. This notifies you that the recommended usage time has been reached. Subsequently the temporary stop notification function does not activate until the main unit power is turned OFF. Depending on the water pressure setting, the level of water remaining in the reservoir after 1 minute may vary.

You can follow the procedure below to turn the notification function ON/OFF.

● Notification function OFF

1. When the main unit power is OFF, press the mode change button once. The main unit starts and the floss mode LED illuminates.
2. Press and hold the mode change button for approximately 3 seconds until the **massage mode** light begins flashing.

● Notification function ON

1. When the main unit power is OFF, press the mode change button once. The main unit starts and the floss mode LED illuminates.
2. Press and hold the mode change button for approximately 3 seconds until the **floss mode** light flashes once.

When the main unit power plug is unplugged, the setting resets to the default setting.

Massage mode

This mode uses a stream of water to stimulate the gums. It is recommended that you follow the procedure below and use this mode in combination with floss mode.

After using floss mode for around 1 minute, press the mode change button and change to massage mode.

Use the water pressure control dial to adjust the water stream to the optimal strength.

When the main unit power is OFF, you can check the current mode from the LED mode indicators by pressing the mode button once. To change the mode, press the mode change button again. You can also change the mode while the product is in use by pressing the mode change button.

When Finished

Turn the power OFF. Empty any liquid left in the reservoir. Unplug the unit.

Caution

When using mouthwash or a solution other than water, read the "Using mouthwash or other solution (P.11)" section.

USING THE TIPS

PIK POCKET™ TIP

The Pik Pocket™ tip is specially designed to direct a jet of water or solution deep into the pocket around the tooth.

When using the Pik Pocket™ tip, **set the main unit water pressure to the lowest setting** and angle the tip at 45 degrees from the tooth. Slowly insert the tip into the pocket around the tooth.

Tongue cleaner tip

Set the water pressure to the lowest setting and place the tip near the center of your tongue. Gently press the tip against your tongue and pull it toward you. Adjust the water pressure to suit your preferences.

Orthodontic tip

When using the orthodontic tip, move the tip slowly and smoothly along your gums, then stop it for a moment to gently clean the space between the teeth and around the orthodontic bracket. Then move on to the next tooth.

Toothbrush tip

Insert the toothbrush tip into your mouth and place the brush head against your gums. It is not necessary to put toothpaste onto the toothbrush tip. Use **a gentle force** and move the brush back and forth. (Use a gentle force that does not bend the bristles, and use short strokes like a manual toothbrush.)

Plaque Seeker™ tip

When using the Plaque Seeker™ tip, position the tip close to your teeth so that the bristles gently contact the teeth.

Move the tip slowly and smoothly along your gums, stopping it briefly at the spaces between the teeth. Brush gently and direct the water stream between the teeth.

Using mouthwash or other solution

The Waterpik® Water Flosser can use mouthwash or other antibacterial solution. In order to prevent clogging, after using a solution other than water, fill the reservoir halfway with water or warm water, then direct the tip toward the sink and operate the main unit until the reservoir is empty.


Storing the tips


There is a tip case which can hold 2 tips located under the reservoir lid. You can easily add or remove a tip by removing the reservoir.

Caution

Not all of these tips are provided with all models.

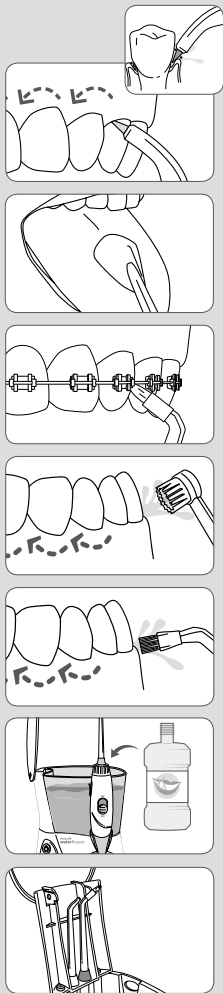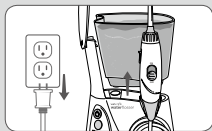Caring for the product (daily)

● Main unit

Turn the power OFF and unplug the power cord from the outlet. Use a soft cloth and a neutral detergent that does not contain a polishing agent, and wipe the product as necessary.

● Reservoir

Pour out the water in the reservoir and dry it in a well ventilated location. If there is dirt on the surface, wipe with a soft cloth and a neutral detergent that does not contain a polishing agent.

The reservoir can be washed in a dishwasher if it is placed on the top rack. (The tips cannot be washed in a dishwasher.)

● Tips

Wipe with a soft cloth and a neutral detergent that does not contain a polishing agent.

Do not wipe the main unit with alcohol, paint thinner, benzene, nail polish remover or other solvent. These may damage the product surface.

- Do not immerse the main unit in water.

Cleaning inside the main unit (once every 1 - 3 months)

Depending on the mineral content of the water that is used, water scale may build up inside the main unit.

If it is not removed, it may decrease product performance.

1. Fill the reservoir with water or warm water and add 2 tablespoons of white vinegar or grain vinegar.
2. Direct the handle and tip toward the sink and operate the main unit until the reservoir is empty.
3. Fill the reservoir with clean water or warm water and repeat the same procedure around 3 times to rinse out the main unit.

Storing the product

- After use, be sure to fully wipe off any water and dry the product before storing it. If the product is left with water on it, it may cause corrosion inside the main unit.
- If the product will not be used for a long period of time, clean the main unit before storing it.
- If the product will not be used, unplug the power plug.
- Do not place the product close to fire, exposed to direct sunlight, in a bathroom or other location exposed to high temperatures and humidity or in an unstable location. Store it in a clean location with good ventilation.
- Store the product in a location out of the reach of children.

Maintenance and inspection

- If the product was not used for a long period of time, be sure to check that the product operates correctly and safely before using it.
- If there is any problem with operation of the product, before requesting repairs, refer to the "Troubleshooting" (P.14) section and check the trouble. If the problem is still not corrected, contact the retailer where you purchased the product, or prepare the purchase invoice and contact the YA-MAN Call Center (P.15).
- If the product has malfunctioned or is damaged, contact the retailer where you purchased the product.

Water leakage

If water leakage occurs as a result of damage to the hose, contact the YA-MAN Call Center (P.15).

TROUBLE SHOOTING / FAQ

Symptom	Possible Causes	Solutions
Unit won't turn on?	The power plug is not inserted into the outlet.	<ul style="list-style-type: none"> - Insert the power plug into the wall outlet.
Low or No water pressure	The tip is clogged, intake port (rubber packing that connects the reservoir and main unit) is clogged, an additive was used, cleaning was insufficient or preparation for use was insufficient.	<ul style="list-style-type: none"> - Perform cleaning. - Check that the reservoir is correctly installed onto the main unit. Gently tilt the reservoir back and forth, securely lower the reservoir lid, and turn ON the main unit power. - Check that the water stream ON/OFF switch is not stuck.
Tip Broken in Unit	The product was dropped.	<ul style="list-style-type: none"> - Screw a metal screw into the end of the broken tip. - Hit eject button and remove screw. Tip will come out with the screw.
Unit is leaking from the hose	A solution other than a designated solution was used (corrosion caused by chemicals), mineral scale accumulated due to insufficient cleaning, or the hose deteriorated when the product was not used for a long time.	<ul style="list-style-type: none"> - Perform cleaning or contact the YA-MAN Call Center (P.15).
Unit is leaking from the base	The product was left with water in it, or the hose is not connected correctly to the main unit.	<ul style="list-style-type: none"> - Drain the water from the reservoir. - Check whether the threaded part where the hose is attached is fully tightened.
Tip Ejects while in use	There is dirt on the tip, or the tip or main unit has malfunctioned.	<ul style="list-style-type: none"> - Try a different tip to see if problem still occurs - If you used additives like mouthwash or peroxide, you will want to flush water in the hole where tip is inserted.

FAQ

Where can I buy replacement tips and accessories? Contact the YA-MAN Call Center (P.15).	How long does the water last? How much does the reservoir hold? The reservoir capacity is approximately 651 mL (enough for approximately 90 seconds of use).
What is the warranty on this product? 2 Year	What is the pressure range? 10-100 PSI, 10 pressure settings
Can I use solutions in my water flosser? If additives (such as mouthwash) are used, rinse the system with a reservoir full of clean water after each use. Failure to rinse the unit may result in loss of product performance and/or function. Growth of bacteria could also be enhanced. Note: iodine, saline or water insoluble concentrated essential oils cannot be used.	My buttons are stuck on the handle what do I do? Soak your handle in a solution of 50/50 water and vinegar for at least 30 minutes. Most likely caused by not cleaning with additive use.
	What is a guideline for how often I should use the product? Up to 2 times/day. You can use it every day.

LIMITED TWO-YEAR WARRANTY

Disposing

When disposing of this product, please dispose of it according to the rules of your area.

Contact Information

If you have any questions about the product warranty, repairs or other matter, please contact the Call Center below.

YA-MAN Call Center

Hours: 9:00 - 17:00 (closed on Saturdays, Sundays, holidays, and during the year-end and New Year holidays)

from standalone or similar phone

 **0120-776-282** (Call is free.)

from a mobile phone or PHS

 **0570-550-637** (A call fee is charged.)

* When requesting repairs, prepare the warranty certificate or purchase invoice that contains the date of purchase and retailer name.

* The Call Center lines are often busy on Mondays.

Product specifications

Rated power	AC 100-240 V, 50/60 Hz
Power consumption	18 W
External dimensions	Approximately W141 x D117 x H263 (mm) when tip is installed
Product weight	Approximately 662 g
Reservoir capacity	Approximately 651 mL (enough for approximately 90 seconds of use)
Accessories	Jet tips: 3, Plaque Seeker tip: 1, orthodontic tip: 1, Pik Pocket tip: 1, Tip Case: 1, instruction manual (with warranty certificate)
Material	Main unit: ABS plastic, reservoir: copolyester, tips: copolyester/ABS plastic
Country of manufacture	China

Warranty information

A warranty for 2 year from the date of purchase is provided and covers problems that occur during ordinary use. (It does not cover consumable supplies or accessories.)

- After the warranty period expires, product shipping expenses will be borne by the customer.
- In the case of problems occurring after the warranty period expires, the customer will be charged for the cost (actual expense) of repairs.
- There are cases when a charge for repairs will be incurred even during the warranty period. Please read the explanation of the warranty terms carefully.

- In the event of a problem with the product, contact the retailer or the YA-MAN Call Center (listed above).

Warranty Terms

1. Problems which occur under conditions of normal use as prescribed in the instruction manual, labels attached to the product, and other precautions shall be repaired at no charge by the retailer where the product was purchased.
2. To receive free repairs of a problem that occurred during the warranty period, present the product and this warranty certificate, and request repairs by the retailer where the product was purchased.
3. If your address will be changed, please consult in advance with the retailer where the product was purchased.
4. If the product was a gift and it is not possible to request repairs by the retailer where the product was purchased listed in this certificate, please contact the YA-MAN Call Center.
5. Even during the warranty period, a charge shall be required for repairs in the following circumstances.
 - (a) The problem or damage is the result of incorrect use or unauthorized repair or modification.

* This certificate is a pledge to provide free repairs based on the warranty period and terms printed on it. Therefore this warranty certificate does not limit the rights of the customer under the law. If you have any questions regarding repairs after the end of the warranty period or other matters, please contact the retailer where the product was purchased.

* A charge will be incurred for repairs in cases of deterioration or consumption resulting from use, or for replacement of consumable supplies.

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**Water Pik, Inc. a subsidiary